

Relevo™

Relevo is a sterile, conformable, superabsorbent dressing for use on exuding wounds. Relevo absorbs and retains exudate in the dressing. The exudate is absorbed vertically by the dressing and retained in the core forming a gel, thus minimizing the risk of maceration and maintaining a moist wound environment.

Relevo


The design allows the dressing to be applied with either side facing the wound. This feature also makes it possible to stack dressings on top of each other if extra absorption capacity is required.

The dressing remains dry and comfortable for the patient. Simultaneously it contributes to a moist wound environment. The hydrophilic contact layer transports fluid effectively away from the wound and distributes it evenly into the ABSORBEST core.


The ultra-thin and highly conformable features allows the dressing to be used on all parts of the body. Relevo is designed to work very well together with compression therapy.




- Cost effective
- User friendly
- Increased patient comfort
- Improved wound healing

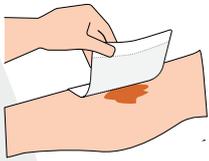


The dressing must be placed close to the wound, or placed in close combination with a suitable dressing.

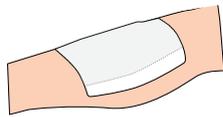

Secure the dressing with suitable bandage, film or tape.

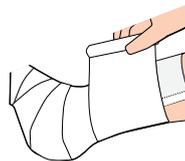

Relevo can be used in conjunction with compression therapy.

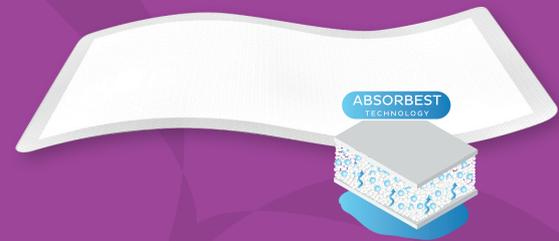

Item	REF	Pcs/box	Product size (cm)	Absorption capacity*
Relevo 10x10 cm	F50011/10	10	10x10	50 ml
Relevo 10x10 cm	F50011/50	50	10x10	50 ml
Relevo 10x20 cm	F50021/10	10	10x20	130 ml
Relevo 10x20 cm	F50021/50	50	10x20	130 ml
Relevo 20x20 cm	F50041/10	10	20x20	250 ml
Relevo 20x20 cm	F50041/50	50	20x20	250 ml



CE 0413

STERILE EO



*) According to international standard EN 13726-1

ABSORBEST

Manufactured by Absorbest AB Sweden
www.absorbest.se | info@absorbest.se