

DONJOY[®]

Performance & Lifestyle

OA Reaction Web[™]



Don't let OA knee pain
slow your patients down!

 **DJO**[®]
GLOBAL

OA Reaction™ Web

A Discreet Solution for Osteoarthritis Knee Pain

- OA Reaction Web™ gives those with OA knee pain the relief they need in a **comfortable, lightweight, discreet solution**
- OA Reaction Web™ knee brace is a **unique and effective way to treat osteoarthritis pain**, featuring patented silicone web technology and comfortable off-loading


Indications

- Early to mild unicompartmental osteoarthritis of the knee & patella
- OA Reaction Web™ is the ideal “first-time” choice: it allows patients to get accustomed to wearing a brace, making it easier to transition into a more rigid OA brace when the time comes
- OA Reaction Web™ is also a good option for patients wishing to use it as a lighter alternative to their usual rigid brace

Patient profile

- Mild to moderate activity levels
- All age groups
- Overweight patients

The technology behind the pain relief

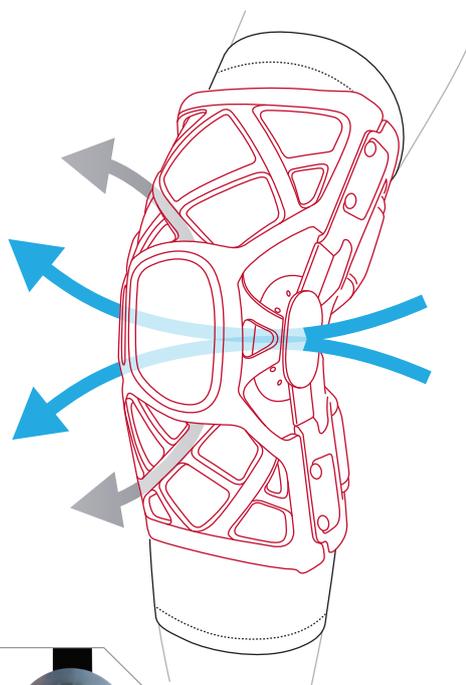

Single malleable aluminium upright

- Built with 6 degrees of valgus correction
- Can be bent to adjust the amount of off-load, to help relieve pressure on painful knee compartment

Off-loading mechanism

- The combined actions of the valgus upright and the condyle pad help mechanically shift the compressive load away from the affected compartment
- Shifting the load to the healthy compartment helps alleviate OA pain

Patented silicone “Web” technology

- Medical grade silicone extends and retracts with leg movement, absorbing and releasing kinetic energy
- Provides shock absorption, anterior knee pain relief and mild medial compartment off-load

Patella stability

- The shape of the anterior cut-out is specifically designed to help realign the patella in the femoral groove and reduce discomfort

Patented, integrated silicone condyle pad

- Provides gentle, comfortable off-loading of the affected compartment
- The hexagonal design on the condyle pad creates an air pocket interface for a non-slip effect



4-Way stretch wraparound design

- Comfortable, form-fitting, easy to apply
- Provides a custom-like fit for most leg shapes

Back pad & straps

- Contoured design & moisture wicking material

Silicone condyle pad

- Gentle off-loading of the knee with a nonslip effect

Padded cuff liners

- For protection and comfort

Single aluminium upright

- Provides stability and can be adjusted to increase off-load

Silicone Web technology

- Provides shock absorption and return of kinetic energy
- The intimate fit increases proprioception, helping patients regain their confidence and stay active

“Open” framework

- The carefully positioned silicone membranes provide patella stability and mild off-load

Low profile, lightweight, discreet

- The level of comfort is such, that patients may forget they are wearing it
- The brace can be easily worn under clothing before, during and after activity

Fitted undersleeve

- For enhanced comfort



OA Reaction™ Web

ORDERING INFORMATION

| Part Numbers | Description |
|--------------|---|
| 82-7426-x | OA Reaction Web - Right Medial/Left Lateral |
| 82-7427-x | OA Reaction Web - Left Medial/Right Lateral |
| 25-9504-x | OA Reaction Softgoods KitRight Med/Left Lat |
| 25-9505-x | OA Reaction Softgoods KitLeft Med/Right Lat |

SIZING INFORMATION

| Brace Size | Thigh Circumference 6" Above Mid-Patella | |
|--------------|--|--------------|
| x = 1 = XS | 33 – 39 cm | 13 – 15.5" |
| x = 2 = S | 39 – 47 cm | 15.5 – 18.5" |
| x = 3 = M | 47 – 53 cm | 18.5 – 21" |
| x = 4 = L | 53 – 60 cm | 21 – 23.5" |
| x = 5 = XL | 60 – 67 cm | 23.5 – 26.5" |
| x = 6 = XXL | 67 – 75 cm | 26.5 – 29.5" |
| x = 7 = XXXL | 75 – 81 cm | 29.5 – 32" |

At DJO Global, we are focused on providing tailored solutions for today's and tomorrow's healthcare challenges, like Osteoarthritis (OA). We believe **Motion is Medicine™** helping patients feel better and regain a healthier lifestyle.

The OA Reaction Web from DonJoy is a new discreet solution for relieving the pain associated with knee OA, while allowing patients to regain or continue their activities.

DJO Global is committed to developing innovative products which are ideally suited to the patients' activity levels.

Discover all our OA solutions at
www.DJOGlobal.com/Motion-is-Medicine

Learn more about OA REACTION WEB at
www.DJOGlobal.com/OAreactionknee

DJO Global | 1a Guildford Business Park | Guildford | Surrey | GU2 8XG | ENGLAND

www.DJOGlobal.eu



Together in Motion™