

Viking™ L mobile lift

For additional strength



With a lifting capacity of 250 kg, Viking L is ideal for most common lifting needs, and offers a flexible lifting solution. The extensive application area means that for lifting medium-weight or larger patients, the Viking L is suitable in most situations, e.g. lifting from the floor, to/from the toilet, horizontal lifting, and/or gait training.

Features and benefits

Optional quick release hook system means caregivers can easily and rapidly change a sling bar, add a Liko Stretch for horizontal lifting or install a scale

Ergonomic hand control for lift operation displays both battery usage and overload indicators

Two lifting speeds accommodate a variety of different lifting situations

User display highlights device usage, including lift count, service required indicator and overloads


Optional armrests help **support patients during ambulation** as well as making it easier for caregivers to manoeuvre the lift

An extensive assortment of accessories is available including slings, slingbars and scales

Viking L mobile lift



Handles

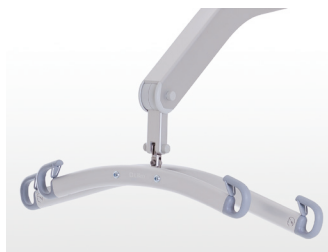

Optional Universal TwinBar 670 slingbar



Hand control



Control box

Technical Specifications

Maximum Load	250 kg
Weight	36.7 kg
Height	1390 to 2030 mm
Protection class	IPX4
Slingbar options	
standard	Universal Slingbar 450
optional	Universal TwinBar 670, Crossbars
Base	
height	140 mm
length	1360 mm
width	695 to 1070 mm
expansion	electric
Castors	
front	100 mm double band
rear	125 mm double band
Material	aluminium
Emergency lifting	electric
Emergency lowering	manual and electric
Lifting speed (without load)	17 mm / sec or 23 mm / sec
Lifting interval	1330 mm
Maximum lifting height	1790 mm

Sound level	51 dB(A)
Battery	lead acid gel or lithium ion
Charger	built in
Hand control	wired
Protection class hand control	IPX4
Product number	2040044

Hill-Rom is a leading global medical technology company with more than 10,000 employees worldwide. We partner with healthcare providers in more than 100 countries, across all care settings, by focusing on patient care solutions that improve clinical and economic outcomes in five core areas: **Advancing Mobility, Wound Care and Prevention, Patient Monitoring and Diagnostics, Surgical Safety and Efficiency and Respiratory Health.** The people, products and programs of Hill-Rom work towards one mission: **Every day, around the world, we enhance outcomes for patients and their caregivers.**

Not all products/options are available in all countries.
For further information about our products and services, please contact your local Hill-Rom representative or visit our webpage:

www.hill-rom.com

Viking L is a versatile lift intended mainly for use in health care, intensive care and rehabilitation.
Class I
Manufacturer: Liko AB, Nedre vägen 100, SE-975 92 Luleå, Sweden

This medical device is a regulated health product which, pursuant to such regulation bears a CE mark. Hill-Rom recommends that you carefully read the detailed instructions for safe and proper use included in the documents accompanying the medical devices. The personnel of healthcare establishments are responsible for the proper use and maintenance of these medical devices.

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

©2017 Hill-Rom Services, Inc. ALL RIGHTS RESERVED.
Doc. No: 5EN183308-02, 12 May 2017

Enhancing outcomes for
patients and their caregivers®

Hill-Rom