



MIDWIFERY

CATALOGUE | 2019



THERE ARE FOUR PARTS TO OUR BUSINESS.

CUSTOMERS | SUPPLIERS | COLLEAGUES | COMMUNITY

ALL OF WHICH CONTRIBUTE EQUALLY TO
OUR ABILITY TO DELIVER AND OUR OVERALL
SUCCESS.

SEAMLESS SUPPORT FOR OUR CUSTOMERS.

Since 1984, USL has grown to be a significant supplier in the healthcare market of New Zealand. We supply a wide range of products to professionals nationwide. We are suppliers to all 20 district health boards and public hospitals in New Zealand. Today we have a significant presence in both consumer retail and industrial markets.

USL Aesthetics is the newest addition to the USL family, with an extensive line of products for the aesthetically minded. By keeping with the USL philosophy we are known for, we are excited about this new direction and look forward to cultivating new partnerships and success in this market.

OUR SERVICING DIVISION.

Our highly skilled service technicians are equipped to install, service and repair the full USL product range and other medical equipment. This includes lasers as well as other servicing and preventative maintenance requirements. USL technicians work according to AS/NZS 3551:2012 standard and our workshop complies with ISO9001.

We are also the only authorised service agent for warranty repairs for all Welch Allyn® products.

IT'S ALL ABOUT PEOPLE.

USL has been built by creating seamless relationships. Since 1984, we've built enduring partnerships with customers, suppliers and our community.

We set out to create the very best environment for our team so that they might deliver their unique skills to our markets.



TABLE OF CONTENTS

HUNTLEIGH SONICAID DIGITAL/FETAL DOPPLER.....	4 - 5
HUNTLEIGH BD4000 AND FM830 CTG MONITORS.....	6 - 7
HUNTLEIGH SONICAID TEAM 3.....	8
HUNTLEIGH FREEDOM.....	9
HUNTLEIGH PULSE OXIMETER.....	10
WELCH ALLYN® MANUAL BP DEVICES.....	11
WELCH ALLYN® LED OPHTHALMOSCOPE & THERMOSCAN.....	12
HILL-ROM® AFFINITY® BIRTHING BED.....	13
PRICE INVENA® NEW PRODUCT.....	14 - 15
USL SPORT SOFT GOODS.....	16 - 17
BIOSKIN™ BRACES AND SUPPORTS.....	18 - 19
BILISCAN & SECA BABY SCALE 354.....	20
CONTINENCE PADS AND UNDERPADS.....	21
SIGVARIS COMPRESSION STOCKINGS.....	22
SIGVARIS VARICOSE VEINS.....	23
SIGVARIS MORNING SICKNESS.....	24
MAM MANUAL BREAST PUMP.....	25
SUDOCREM® HEALING CREAM.....	26
INFACOL® AND SEABAND® MAMA!.....	27
MILTON® AND D-BRIEF.....	28
PRODUCTS AND PRICING.....	29 - 31
CONTACT DETAILS.....	32



HUNTLEIGH SONICAID DOPPLER

THE SONICAID DIGITAL DOPPLER

A NEW DIMENSION IN SOUND AND VISION

Visual presentation of fetal heart rate (FHR), in numeric and trace views, providing added clinical reassurance.

With over 40 years' experience in fetal monitoring, our technological innovations have enabled us to develop a new digital platform for our latest doppler range.

This generation of doppler combines excellence in probe performance with digital signal processing, to provide superb digital audio using our unique Dynamic Digital Noise Reduction (DDNR) system and gel filter.

Our latest digital algorithms optimise FHR performance and the high resolution colour display includes new 'BIG numbers' & trace display modes. Doppler sound recording, on-screen trace review and playback are just some of the many features in this state-of-the-art digital doppler unit.



\$1,199.00 + GST
Normally \$1,610.00 + GST
Code: 14646

HUNTLEIGH FETAL DOPPLERS

FD1 FETAL DOPPLER

The FD1 combines all the powerful features of the D920 with the added advantage of displaying fetal heart rate. The high sensitivity 2Mhz probe provides reliable detection throughout pregnancy.

- Advanced fetal heart rate (FHR) processing
- Active noise reduction for crystal clear fetal heart sounds
- 2Mhz waterproof probe
- Advanced power management to conserve battery life
- Easy clean probe for improved infection control



D920 FETAL DOPPLER

This audio doppler has been designed with ease-of-use in mind. Its high sensitivity transducer coupled with clear fetal heart sounds, enable the clinician to quickly locate the fetal heart.

- High sensitivity ultrasound for improved ease of use
- 2Mhz waterproof transducer suitable for use in a water birthing environment
- Battery low indicator
- Easy clean transducer for improved infection control



HUNTLEIGH BD4000 CTG MONITOR



ENTRY LEVEL FETAL MONITORING

The BD4000xs series fetal monitors provide solutions for all areas of fetal monitoring. Using state-of-the-art digital signal processing technology with enhanced FHR performance, the range includes antenatal and intrapartum models and a range of options and interface capabilities.

This product provides the most cost effective and flexible solution for all your fetal monitoring needs. Its compact size and styling belie its ruggedness and robust design, using the same materials, construction and design expertise as our world renowned handheld dopplers.

FEATURES:

- Twins capability as standard
- Unique no compromise full size twin FHR scale printing
- User adjustable, tachycardia, bradycardia, loss of contact and cross-channel rate verification alerts

MONITOR INTERFACE:

- User-customisable trace annotation facility
- Automatic fetal movement detection
- Intuitive user interface and controls simplifying training
- Maternally sensed fetal movement marker
- FECG and IUP upgrade options



P.O.A

Please contact us on
0800 658 814 for further
pricing



HUNTLEIGH
FM830E CTG MONITOR

HIGH RISK FETAL MONITORING

"Great performance fetal monitoring even when used on women with a high body mass index (BMI)...Reduces the need for invasive procedures"

CLINICIANS REPORT A WIDE RANGE OF BENEFITS INCLUDING:

- Improved monitoring experience for mothers
- Reduced need for invasive FECG monitoring
- Reduced costs
- Reduced risk of infection and cross contamination

FEATURES:

- Maternal vital signs
- Quick view trace display
- High sensitivity transducers
- Unique fetal CTG analysis
- Unique TREND intrapartum recording
- Ergonomic design

HIGH SENSITIVITY TRANSDUCERS

New advances in biocompatible materials have enabled us to develop the most sensitive uterine activity monitoring for patients of all sizes.

Waterproof design makes for thorough and easy cleaning; essential for infection control.



P.O.A

Please contact us on
0800 658 814 for further
pricing

SONICAID TEAM 3

THE NEW SONICAID TEAM 3

FETAL MONITORING FOR THE NEXT GENERATION

Stylishly designed not to dominate the bedside, Team 3 offers accurate and clear visual status of the fetus and mother. The range is full of features providing cost-effective, reliable and accurate fetal/maternal monitoring.

Easy to use and intuitive via the icon driven touchscreen, the fetal heart rate can be displayed as 'Big numbers' which auto scale for single, twins or triplet monitoring for optimum visibility.

With the widest range of options, Team3 provides excellence in fetal heart rate monitoring with high sensitivity transducers and unique digital processing delivering the best possible monitoring of all mums, from low BMI to high in all antenatal and labour settings.

P.O.A

Please contact us
on 0800 658 814 for
further pricing



TEAM3 - A/P



TEAM3 - I/P

Antepartum and Intrapartum models available in singleton, twins or triplet configuration. Available with maternal vital signs and integral battery options.



TEAM3 - E-CTG

Paperless monitoring option with no printer.

HUNTLEIGH



HUNTLEIGH
FREEDOM

WIRELESS FETAL MONITORING SOLUTION

Sonicaid Freedom transducers offer flexibility of choice for monitoring mother and baby. The system consists of wireless fetal heart rate and contractions transducers and a receiver base system.

Sonicaid Freedom is designed for use with the FM800 Encore series of fetal monitors from Huntleigh.

Our waterproof wireless transducers allow the fetus to be continuously monitored while mum is bathing, showering or on the move, and provides the widest choice of birthing positions.



P.O.A

Please contact us
on 0800 658 814 for
further pricing





HUNTLEIGH PULSE OXIMETER

MINIPULSE HANDHELD PULSE OXIMETER

The professionals choice, ideal for both spot checking and continuous monitoring of adult, paediatric and neonatal patients.

The MiniPulse range has been designed to provide clinicians with fast, reliable, portable measurements in almost any setting.

FEATURES:


- Larger display
- Long battery life
- Value for money
- Full function alarms
- Robust casing
- Probe storage ability
- Easy to use




INFANT WRAP
(3-15kg)

\$250.00 + GST

Code: 10383



**PAEDIATRIC
PROBE (5-15kg)**

\$240.00 + GST

Code: 10387

MINI PULSE
\$950.00 + GST

Normally \$1,127.50 + GST
Includes adult probe
Code: 10366



WELCH ALLYN®
MANUAL BP DEVICES

WELCH ALLYN® MANUAL BP DEVICES

- Dura shock technology
- Remain in calibration longer than traditional gauges so is more accurate
- The only gauge that passes the AAMI shock resistant standard - can be dropped 76cm and remain in calibration
- Latex free
- Warranty - 3 years on cuffs

OPTION TWO

DS55 has a traditional thumbscrew-style air release valve - 5yr warranty on calibration

\$170.00 + GST

Code: 5041

OPTION ONE

S66 manometer and adult cuff in zip up case - 15yr warranty on calibration

\$191.00 + GST

Normally \$250.00 + GST
Code: 40509827





WELCH ALLYN® LED OPHTHALMOSCOPES

- Brand new technology
- Pocket Plus includes case
- Available in colours below, limited numbers at this price
- Warranty – 5 years



LIMITED OFFER
\$240.00 + GST
 Normally \$265.00 + GST

Prices subject to change
 with no notice



CASE INCLUDED

WELCH ALLYN® EAR THERMOMETER

The PRO 6000 advantage lies in our innovative PerfectTemp technology, which overcomes the potential for low readings compared to core, by adjusting for factors that impact measurement accuracy.

Designed to take quick and accurate readings, to improve patient satisfaction.

- Large, backlit LCD display
- Dedicated memory recall button to display last measurement taken
- 60 second timer for manual pulse and respiration measurement
- Pre-heated probe tip and sensor helps support accurate measurements
- Green illuminating LED signals ExacTemp positioning feedback on probe tip
- Ergonomic handle design for a comfortable and secure hold



\$655.00 + GST
 Normally \$310.00 + GST
 Code: 406000200

Prices subject to change
 with no notice



HILL-ROM® AFFINITY® BIRTHING BED

As used by all but one DHB in New Zealand

Features

- Straight edge (U-Cut) mattress
- Integrated calf supports electric variable height
- Head and foot elevation
- Two half length side rails with integrated patient control
- Central brake and steer
- CPR
- Removable headboard
- Battery back-up
- IV Pole
- Automatic night light for patient security
- Trendelenberg manual operation of functions
- SWL 227kg
- Used by 19 DHB's – including Canterbury DHB



P.O.A

Please contact us
on 0800 658 814 for
further pricing

NEW PRODUCT RANGE

Exclusive to **usl** medical

Available in
**APRIL
2018**

CLAMPCUT®

All-in-one umbilical cord cutter



ClampCut all-in-one is a new innovative device for clamping and cutting the umbilical cord during childbirth.

- ✓ Reduces blood spatter
- ✓ Umbilical cord is fixed
- ✓ No chance of sharp injures

CLAMCORD®

Three-legged umbilical cord clamp



ClamCord is an umbilical cord clamp with three legs, preventing the cord from slipping out during use.

- ✓ Umbilical cord is fixed
- ✓ Easy to close
- ✓ Loud “click” when closed

CLAMPOFF®

Umbilical cord clamp remover



ClamOff is an universal umbilical cord clamp remover for all types of umilical cord clamps.

- ✓ Can remove all clamps
- ✓ Easy to close
- ✓ Ergonomic handle

CLAMPCUT®

All-in-one umbilical cord cutter



ClampCut® All-in-one

ClampCut® is a sterile disposable obstetric instrument for clamping and cutting the umbilical cord at child birth.

The intention with ClampCut® is to make births more safe and simple for both the parents, midwives, doctors and nurses.

Once you have experienced ClampCut®, you will never return to the traditional procedure for clamping and cutting the umbilical cord.

Safe

ClampCut® has been designed to minimize blood spatter and the chance of sharps injuries. It is made of nylon and do not contain latex.

Simple

The traditional clamping and cutting procedure with six instruments is replaced with one easy 'Click-Click-Cut' action.

Superior

ClampCut® has a triangular cavity, which secures the umbilical cord before it is clamped and then cut. This stabilizes the clamping and cutting action.

Sustainable

ClampCut® is an environment-friendly device. It's plastic components are manufactured from nylon which evaporates to water when incinerated.

ClampCut® is the safe choice

Easy for the partner to cut the cord

No chance of blood spatter and spreading of blood born diseases

Accountable for 2 units only - one on the baby and one on the placenta

How ClampCut® works



Place the umbilical cord in the clamps. Make sure the white clamp is facing the infant. Close the clamp.

Close the green clamp facing the placenta. Now the umbilical cord is ready to be cut.

Press the handle on the green clamp and push both clamps through the cutting house.



Traditional birth set

CLAMCORD®

Three-legged umbilical cord clamp

ClamCord®

ClamCord® is a sterile disposable obstetric instrument for clamping the umbilical cord at child birth.

The intention with ClamCord® is to make births more safe and simple for both the parents, midwives, doctors and nurses.

Once you have experienced ClamCord®, you will never return to the traditional procedure for clamping the umbilical cord.

Safe

ClamCord® has been designed to minimize wrong placement of cord clamps on umbilical cords. ClamCord® is made of nylon and do not contain latex.

Simple

The traditional umbilical cord clamp has only two legs, which means the cord can slip out during use. The third leg on ClamCord® prevents this.

Superior

ClamCord® is a small umbilical cord clamp and does not annoy the baby during the period where the umbilical cord stump will fall of.

Sustainable

ClamCord® is an environment-friendly device. The plastic body is manufactured from nylon which evaporates to water when incinerated.


ClamCord® is the safe choice

Easy closing of the clamp

Umbilical cord is fixed inside the clamp

Loud audible click when the clamp gets closed

How ClamCord® works

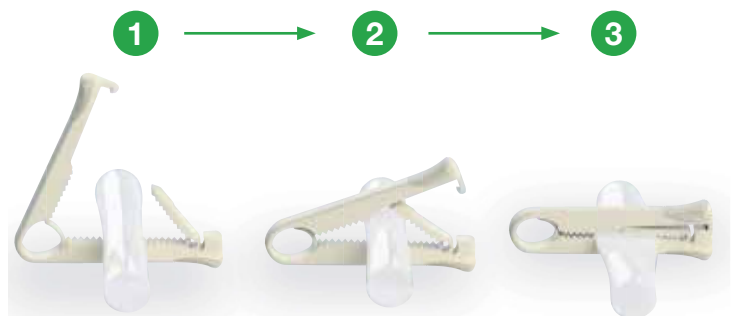

USL MEDICAL BRACES & SUPPORTS

WRIST/PALM SUPPORT

USL Wrist/Palm Support is ideal for post fracture, immobilisation, wrist sprains and strains.

- Breathable polyester construction
- Hook and loop closures
- Malleable aluminium splints
- Available in left or right

WRIST/THUMB SUPPORT

USL Wrist/Thumb Support is ideal for post fracture, thumb immobilisation, wrist sprains and strains and support of the CMC and MCP joint.

- Breathable polyester construction
- Hook and loop closures
- Malleable aluminium splints wrist and thumb
- Available in left or right



CODE	DESCRIPTION	SIZE
LEFT		
10901-L	USL Sport Wrist/Palm Support	X-Small
10217-L	USL Sport Wrist/Palm Support	Small
10218-L	USL Sport Wrist/Palm Support	Medium
10219-L	USL Sport Wrist/Palm Support	Large
10220-L	USL Sport Wrist/Palm Support	X-Large
RIGHT		
10902-R	USL Sport Wrist/Palm Support	X-Small
10217-R	USL Sport Wrist/Palm Support	Small
10218-R	USL Sport Wrist/Palm Support	Medium
10219-R	USL Sport Wrist/Palm Support	Large
10220-R	USL Sport Wrist/Palm Support	X-Large

Code	Description	Size
LEFT		
10221-L	USL Sport Thumb/Wrist Support	Small
10222-L	USL Sport Thumb/Wrist Support	Medium
10223-L	USL Sport Thumb/Wrist Support	Large
10224-L	USL Sport Thumb/Wrist Support	X-Large
RIGHT		
10221-R	USL Sport Thumb/Wrist Support	Small
10222-R	USL Sport Thumb/Wrist Support	Medium
10223-R	USL Sport Thumb/Wrist Support	Large
10224-R	USL Sport Thumb/Wrist Support	X-Large

* See page 15 for USL Wrist Splint size chart

UNIVERSAL THUMB WRAP

USL Sport Universal Thumb Wrap is recommended for partial immobilisation after surgery or trauma.

- Relieves pain and stabilises CMC joint
- Aluminium insert provides thumb stability
- Fits left or right



\$7.80 + GST

CODE	DESCRIPTION	UNIT
10216	USL Sport Universal Thumb Wrap	Each

WRIST/THUMB WRAP NEOPRENE

USL Wrist/Thumb Wrap is an adjustable support for the wrist and thumb.

- Wrap-around design
- Neoprene promotes warmth and healing
- One size fits all
- Available in left or right



\$7.25 + GST

CODE	DESCRIPTION	SIZE
6972	Wrist/Thumb Wrap Neoprene Universal	Left
6971	Wrist/Thumb Wrap Neoprene Universal	Right

SIJ BELT

- Relieves stress and instability
- Can be used during pregnancy or post-partum



\$9.90 + GST

CODE	DESCRIPTION	SIZE
13903	USL SIJ Belt	S/M
13904	USL SIJ Belt	L/XL

WRIST/THUMB STABILISER

USL Wrist/Thumb Stabiliser is recommended for partial immobilisation of the thumb. The aluminium metal stay provides thumb stability.

- CMC joint injury
- Thumb strain/sprain
- Arthritic pain
- De Quervain's Disease

\$7.80 + GST



CODE	DESCRIPTION	SIZE	MEASUREMENT (TAKEN AROUND THE WRIST)
13905	USL Wrist/Thumb Stabiliser	Small	15.25 cm - 19 cm
13906	USL Wrist/Thumb Stabiliser	Medium	17.75 cm - 21.5 cm
13907	USL Wrist/Thumb Stabiliser	Large	19 cm - 23 cm
13908	USL Wrist/Thumb Stabiliser	X-Large	21.5 cm - 25.5 cm
13909	USL Wrist/Thumb Stabiliser	XX-Large	23 cm - 26.75 cm

USL WRIST SPLINT SIZING CHART

MEASUREMENT IS TAKEN AROUND THE WRIST



X-SMALL	10CM-12.5CM
SMALL	12.5CM-15CM
MEDIUM	15CM-18.75CM
LARGE	18.75CM-22.5CM
X-LARGE	22.5CM-26.25CM

BIOSKIN™ BRACES & SUPPORTS

DP2™ & DP3™ COCK-UP WRIST

FEATURES

- Low profile and comfortable
- Soft, fleece material against skin
- Allows full use of fingers and thumbs
- Removable/bendable dorsal and palmar stays
- Available left or right

INDICATIONS

- Carpal Tunnel Syndrome
- Arthritis
- Wrist sprain
- Bursitis



DP2
\$34.50 + GST
DP3
\$40.50 + GST

CODE	DESCRIPTION
BS5423*L	DP2 Cock-Up Wrist Splint 6 1/2" - Left
BS5423*R	DP2 Cock-Up Wrist Splint 6 1/2" - Right
BS5424*L	DP3 Cock-Up Wrist Splint 8 1/2" - Left
BS5424*R	DP3 Cock-Up Wrist Splint 8 1/2" - Right

BOOMERANG

FEATURES

- Universal size
- Easy attachment
- Soft, fleece material against skin
- Fits left or right

INDICATIONS



- General wrist support
- Carpal Tunnel Syndrome
- Arthritis
- Wrist sprain
- Wrist pain
- Bursitis



BOOMERANG
\$32.50 + GST

CODE	DESCRIPTION
BS5420*L	Boomerang Wrist Wrap - Left
BS5420*R	Boomerang Wrist Wrap - Right

WRIST

	CIRCUMFERENCE OF EXTENDED WRIST			
	SIZE	*PN	INCHES	CENTIMETERS
	XS - S	1	4 - 6	10 - 15
	M - L	3	6 - 8	15 - 20
	XL - 2XL	5	8 - 10	20 - 25
	THUMB SPICA			
	Length of thumb			
	SIZE	*PN	INCHES	CENTIMETERS
	XS - M	2	3.25 - 4.25	8 - 11
	L - 2XL	5	4.25 - 5.25	11 - 14

SIZE INDICATION



To indicate size, replace * with corresponding *PN (Part Number) from the sizing chart.

THUMB SPICA

FEATURES

- Controls motion of the CMC and MP joints
- Low profile
- Soft, fleece material against skin
- Fits left or right

INDICATIONS

- De Quervains
- Game keeper's thumb
- Bursitis
- Sprains and strains

THUMB SPICA
\$35.00 + GST



CODE	DESCRIPTION
BS5450*	Thumb Spica

*Refer to size chart on page 18

WRIST/THUMB SPICA

FEATURES

- Controls motion of the CMC and MP joints
- Cock-up wrist support
- Soft, fleece material against skin
- Available left or right

INDICATIONS

- De Quervains
- Scaphoid Fracture
- Bursitis
- Sprains and strains

WRIST/
THUMB SPICA
\$42.00 + GST



CODE	DESCRIPTION
BS5452*L	Wrist/Thumb Spica - Left
BS5452*R	Wrist/Thumb Spica - Right

*Refer to size chart on page 18



STANDARD
ELBOW SKIN™
\$33.50 + GST

STANDARD ELBOW SKIN™

FEATURES

- Full elbow coverage
- Comfortable, low-profile fit

INDICATIONS

- Bursitis
- Joint swelling
- Biceps/triceps strains

CODE	DESCRIPTION
BS4210*	Standard Elbow Skin - Ultima Material



TENNIS
ELBOW SKIN™
\$35.00 + GST

TENNIS ELBOW SKIN™

FEATURES

- Silicone pad fits over painful area
- Full elbow coverage

INDICATIONS

- Medial Epicondylitis (Golfer's Elbow)
- Lateral Epicondylitis (Tennis Elbow)

CODE	DESCRIPTION
BS4220*	Tennis Elbow Skin - Ultima Material



TENNIS
ELBOW BAND™
\$34.00 + GST

TENNIS ELBOW BAND™

FEATURES

- Silicone pad fits over painful area
- Fits comfortably below elbow

INDICATIONS

- Medial Epicondylitis (golfer's elbow)
- Lateral Epicondylitis (tennis elbow)

CODE	DESCRIPTION
BS4250*	Tennis Elbow Band - Ultima Material
BS5250*	Tennis Elbow Band - 2s Material

ELBOW

CIRCUMFERENCE OF LARGEST PART OF FOREARM				
SIZE	*PN	INCHES	CENTIMETERS	
XS	0	8 - 9	20 - 23	
S	1	9 - 10	23 - 26	
M	2	10 - 11	26 - 28	
L	3	11 - 12	28 - 31	
XL	4	12 - 13	31 - 33	
2XL	5	13 - 14	33 - 36	

SIZE INDICATION



To indicate size, replace *with corresponding *PN (Part Number) from the sizing chart.

BILISCAN & BABY SCALES

\$6,600.00 + GST

Normally \$7,659.00 + GST
Code: BSCANHUL

BILISCAN

BiliScan is a portable, non-invasive, transcutaneous bilirubin-analysis instrument.

YOUR ADVANTAGES

- Better measurement method because of innovative technology
- High measuring accuracy because of bichromatic photometry
- High correlation of measured value with plasmabilirubin concentration
- Fast decision - making because of immediate results
- Bloodless, painless and transcutaneous measurement
- Easy handling due to ergonomic design and well thought-out technology, compact and easy to use
- Energy-efficient because of energy-saving programming
- Cost effective in purchase and maintenance
- No consumables needed



SECA BABY SCALES 354

The seca 354 is an all-around scale which has been developed on the basis of practical experience. It can be used not only as a baby scale but also as a floor scale for small children.

- Babies can be weighed safely and comfortably in the tray of the seca 354
- Baby scale and floor scale for children in one
- Mobile and robust for transport
- Breast-Milk-Intake-Function

*Bag not included

SCALES
\$240.00 + GST
Code: 12997

BAG
\$70.00 + GST
Code: 12993





HARTMANN CONTINENCE

MOLICARE LADY PADS PKT 14

These are especially designed for light to moderate incontinence. The anatomically shaped pad can be fitted securely and discreetly in firm fitting underwear or MoliCare Premium Fixpants.

- Soft breathable material for healthy skin
- Sift top sheet treated with Aloe Vera
- Three layer absorbent core that pulls fluid away from surface so that skin remains dry and healthy
- Breathable textile back sheet



MOLICARE PADS PKT 28

Are especially designed for light to moderate incontinence. The anatomically shaped pad can be fitted securely and discreetly in firm fitting underwear or net pants.

- Absorbent pads with a Dry Plus fluid distribution layer and an odour neutralising super absorbent core
- Adhesive strip to secure in underwear or MoliCare Fixpants®



MOLICARE FIX PANTS

Washable or can be single use used for holding pads in place instead of underwear.



MOLINEA UNDERPADS

Protective waterproof backing and soft top sheet. Elastic foil stamping at edges avoids leakage and prevents underpad from moving.

- Useful to protect bed during labour
- Available in two sizes and by carton or packets





SIGVARIS

SIGVARIS COMPRESSION STOCKINGS

Pregnancy is a very special time; give your legs the attention they deserve. Relieve the feeling of tired, heavy, aching legs. Assist with the reduction of swollen feet and ankles. Reduce the risk of varicose veins and increase your quality of life with the stylish SIGVARIS range.

30% of pregnant women will develop varicose veins in their first pregnancy **and 55% in their second.**

WHAT ARE THE RISKS?

It is important to improve varicose vein circulation problems, as blood supply is vital to supply the body and organs with nutrients, oxygen and remove waste products. Varicose veins while unattractive, can also become inflamed and lead to thrombosis and phlebitis. They can also indicate that the body's major organs are receiving a reduced blood supply.


When a pregnant woman's uterus receives a reduced bloody supply this constitutes a stressful situation for the unborn baby. If pregnant woman's feet are already swollen on rising in the morning, a doctor should be consulted.

HOW DO SIGVARIS MEDICAL GRADE COMPRESSION STOCKINGS WORK?


Preventing or treating varicose veins can be achieved through SIGVARIS medical grade compression stockings. SIGVARIS compression stockings are capable of mechanically counteracting the raise internal pressure with controlled external pressure. This reduces the volume of blood in the feet and leg veins driving it back to the heart and major organs.

Wearing compression stockings during pregnancy prevents blood from pooling in the veins and reduces the increase of an elevated heart rate in both mother and child.

REDUCED BLOOD SUPPLY:

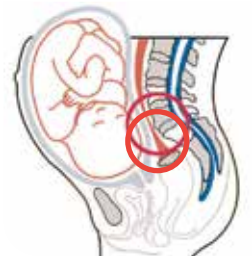


HEALTHY VEIN BLOOD FLOW:

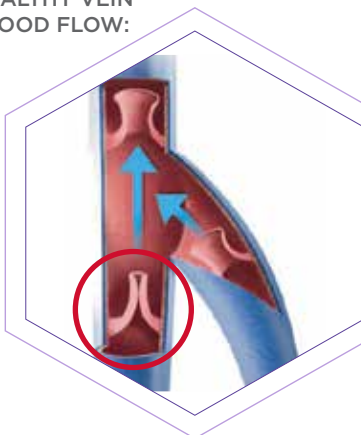

VARICOSE VEINS:



VARICOSE VEIN BLOOD FLOW:

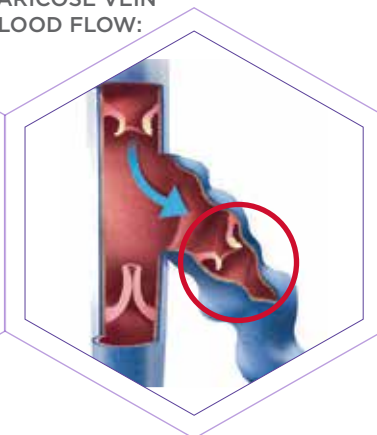

Suffering from morning sickness?

Compression stockings could be a non-pharmaceutical treatment to alleviate nausea and vomiting during early pregnancy.

SIGVARIS

A randomized crossover trial on the effect of compression stockings on nausea and vomiting in early pregnancy



BACKGROUND

50-80% of all women experience nausea and vomiting in early pregnancy (NVP). NVP has a detrimental impact on women's quality of life (QoL), especially their ability to work.

Despite the fact that a number of therapeutic interventions exist to alleviate symptoms – ranging from pharmaceutical products to natural therapies – women remain reluctant to take medications during pregnancy, and a safe, effective and acceptable treatment for NVP is still to be found.

STUDY DESIGN

The aim of this study was to evaluate the effect of medical compression stockings (MCS) on symptoms of nausea, vomiting and dizziness between 8 and 16 weeks of pregnancy, and their impact on QoL factors.

PARTICIPANTS

Pregnant women >18 years and with mild to moderate NVP were eligible to participate. A total of 74 women were enrolled in the study. 58 women completed the study and were included in the analysis set.



GROUP 1: «MCS first» = 2 weeks with // 2 weeks without

- Eligible: n=33
- Dropout: n=4 (3 miscarriages, 1 questionnaire incomplete)
- Final evaluation: n=29

PATIENT RANDOMIZATION



GROUP 2: «MCS second» = 2 weeks without // 2 weeks with

- Eligible: n=41
- Dropout: n=12 (3 miscarriages; 9 did not return for controls)
- Final evaluation: n=29



RESULTS: Endpoint = change in the Nausea and Vomiting in Pregnancy Quality of Life (NVPQOL) score



Compliance rate



Would recommend wearing MCS during pregnancy



CONCLUSION

Wearing compression stockings in early pregnancy improves not only nausea and vomiting-associated symptoms, but also QoL factors such as dizziness, leg pain and physical/physiological limitations.

AUTHORS' SUGGESTION

Compression stockings could be a non-pharmaceutical instrument to alleviate nausea and vomiting symptoms during early pregnancy

* Reference: A randomized crossover trial on the effect of compression stockings in nausea and vomiting in early pregnancy. Mendoza E, Amsler F. Int J Womens Health. 2017;9:89-99.



MAM Breast Care Pump

0+ months 1+ BPA FREE



The Manual Breast Pump from MAM has been developed with real attention to detail and it can be adapted to meet individual needs.

There are three reasons why expressing milk with the Manual Breast Pump is very pleasant and simple. The suction controller allows the suction strength to be adjusted easily to your individual preference.

When it's turned all the way to the right, the vacuum is held for up to three seconds. In this position, lots of milk can be expressed quickly and easily.

Pack includes: 1x 130ml Anti Colic bottle & 2x sealing disks

MAM Breast Care Accessories



Breast Care Pads

Extra thin and very absorbant



Nipple Care Shields

For short term support during breastfeeding



Microwave Steriliser

For microwave and cold water sterilisation



4 Parts
ONLY

Easy cleaning & assembling

Sudocrem®

HEALING CREAM

Happy bum. Happy mum.



Triple Action Emollient Cream Soothes, Heals, Protects.

Sudocrem® helps to soothe, heal and protect your baby's delicate skin and works in three ways:

- A water-repellent base forms a protective barrier, helping to stop any irritants (urine and faeces) coming into contact with the skin
- It both softens and soothes sore inflamed skin
- It helps to ease pain and irritation



Sudocrem®
can be applied
at every nappy
change to help
soothe your
baby's skin.

*Always read the label. Use only as directed. If symptoms persist, see your healthcare professional.



Infacol Wind Drops

RELIEF FROM INFANT WIND

CHOOSE A SILENT NIGHT

Chosen by mums for treating infant wind. Infacol® Wind Drops are an effective method of treating wind in babies and infants, and have been helping parents for generations.

- Infacol has been specially formulated to relieve wind, infant colic and griping pain
- Clinically proven to significantly reduce the frequency and severity of crying attacks associated with colic
- Can be used from birth onwards
- Natural orange flavour
- Sugar, alcohol and colourant free
- Comes with a handy integrated dropper - making it easy to give to your baby
- Each bottle contains enough liquid for up to 100 feeds

SEA-BAND® mama!

Drug Free Nausea Relief Supporting Tender Tummies Through Pregnancy



The new Sea-Band® Mama! range has been specially created to help soothe and calm your tummy during your pregnancy so that you can start enjoying this special time again.

Sea-Bands use acupressure and can be used in conjunction with our Ginger Lozenges to give natural **DRUG FREE** nausea relief.

Drug free, free of artificial colour and flavours, the Sea-Band® Mama! range can support you throughout your pregnancy.



Milton®

TRUSTED SINCE 1916

Cold Water Sterilisation in just 15 minutes!

No boiling, no rinsing, no hassles

Milton® Anti-bacterial Tablets provide hospital grade disinfection in only 15 minutes. Clinically proven to kill 99.9% of germs. No need to boil or rinse and no residual taste remains.

Ideal for disinfecting baby feeding accessories, breast pump parts, bottles, soothers, teething rings, small plastic toys, bowls and plastic cutlery, to name a few.



D»Brief



MENS

WOMENS

CHILDRENS

MEN & WOMEN PADS



Real underwear for light to medium bladder weakness

PRODUCTS

CONSUMABLE MATERNITY

PRICING INDEX

CODE	DESCRIPTION	PRICE EXCL GST	UOM
BRACING & SUPPORTS			
BS54231L	BioSkin DP2 Cock Up Wrist Left	\$34.50	EACH
BS54233L	BioSkin DP2 Cock Up Wrist Left	\$34.50	EACH
BS54235L	BioSkin DP2 Cock Up Wrist Left	\$34.50	EACH
BS54231R	BioSkin DP2 Cock Up Wrist Right	\$34.50	EACH
BS54233R	BioSkin DP2 Cock Up Wrist Right	\$34.50	EACH
BS54235R	BioSkin DP2 Cock Up Wrist Right	\$34.50	EACH
BS54241L	BioSkin DP3 Cock Up Wrist Left	\$40.50	EACH
BS54243L	BioSkin DP3 Cock Up Wrist Left	\$40.50	EACH
BS54245L	BioSkin DP3 Cock Up Wrist Left	\$40.50	EACH
BS54241R	BioSkin DP3 Cock Up Wrist Right	\$40.50	EACH
BS54243R	BioSkin DP3 Cock Up Wrist Right	\$40.50	EACH
BS54245R	BioSkin DP3 Cock Up Wrist Right	\$40.50	EACH
13903	USL Sport Sacroiliac Belt S/M	\$9.90	EACH
13904	USL Sport Sacroiliac Belt L/XL	\$9.90	EACH
10216	USL Sport Universal Thumb Wrap (Elasticated)	\$7.80	EACH
6971	USL Sport Wrist and Thumb Neoprene Right	\$7.25	EACH
6972	USL Sport Wrist and Thumb Neoprene Left	\$7.25	EACH
10217-L	USL Sport Wrist/Palm Support Left Small	\$12.50	EACH
10218-L	USL Sport Wrist/Palm Support Left Medium	\$12.50	EACH
10219-L	USL Sport Wrist/Palm Support Left Large	\$12.50	EACH
10220-L	USL Sport Wrist/Palm Support Left X-Large	\$12.50	EACH
10901-L	USL Sport Wrist/Palm Support Left X-Small	\$12.50	EACH
10217-R	USL Sport Wrist/Palm Support Right Small	\$12.50	EACH
10218-R	USL Sport Wrist/Palm Support Right Medium	\$12.50	EACH
10219-R	USL Sport Wrist/Palm Support Right Large	\$12.50	EACH
10902-R	USL Sport Wrist/Palm Support Right X-Small	\$12.50	EACH
13905	USL Sport Wrist/Thumb Stabiliser Universal S	\$7.80	EACH
13906	USL Sport Wrist/Thumb Stabiliser Universal Medium	\$7.80	EACH
13907	USL Sport Wrist/Thumb Stabiliser Universal Large	\$7.80	EACH
13908	USL Sport Wrist/Thumb Stabiliser Universal X-Large	\$7.80	EACH
13909	USL Sport Wrist/Thumb Stabiliser Universal XXL	\$7.80	EACH

CODE	DESCRIPTION	PRICE EXCL GST	UOM
continued. BRACING & SUPPORTS			
10221-L	USL Sport Wrist/Thumb Support Left Small	\$7.80	EACH
10221-R	USL Sport Wrist/Thumb Support Right Small	\$7.80	EACH
10222-L	USL Sport Wrist/Thumb Support Left Medium	\$7.80	EACH
10222-R	USL Sport Wrist/Thumb Support Right Medium	\$7.80	EACH
10223-L	USL Sport Wrist/Thumb Support Left Large	\$7.80	EACH
10223-R	USL Sport Wrist/Thumb Support Right Large	\$7.80	EACH
10224-L	USL Sport Wrist/Thumb Support Left X-Large	\$7.80	EACH
10224-R	USL Sport Wrist/Thumb Support Right X-Large	\$7.80	EACH
CONSUMABLES			
14215	Cotton Buds Double End Plastic Stick	\$4.60	Box 200
9446	Couch Roll Pure Cellulose 49cm x 50M	\$15.00	Roll
10275	HealthE Hand Sanitiser Gel 375ml	\$8.20	EACH
20001	HealthE Hand Sanitiser Gel 60ml	\$3.20	EACH
8201	Measure Tape Baby Disposable	\$32.95	Pkt 500
12857	Reynard Surface Disinfectant Wipe 33x22cm	\$5.20	Pkt 50
12858	Reynard Surface Disinfectant Wipe Bucket 28x25cm	\$18.85	Pkt225
15703	Reynard Wipes IPA Surface Disinfecting 23x14cm	\$13.76	Pkt 150
4614	Sharps Container BD 1.4 Litre	\$11.53	EACH
14214	Swab Mouth USL Cotton Tip 20cm	\$3.45	Pkt 50
12955	USL Alcohol Swab 70% Isopropyl Alcohol	\$4.50	Box 200
4827A	Lubrigel 500gm Pump Pot	\$7.80	EACH
12859	Reynard Surface Disinfectant Spray 1 Litre	\$20.95	EACH
2387	Sodium Chloride .9% Injection 10ml	\$27.55	Box 50
11581	Surgi-Gel Plus Lubricating Sachets 3ml	\$22.00	Box144
11972	Aquasonic Clear Gel 5L	\$ 68.00	EACH
DIAGNOSTICS			
406000200	Pro 6000 Thermoscan W/A	\$265.00	EACH
10047	Huntleigh Fetal Doppler with LCD Screen - FD1	\$950.00	EACH
10046	Huntleigh Fetal Doppler - D920	\$799.00	EACH
14646	Huntleigh Fetal Doppler SR2 (including probe)	\$1,199.00	EACH
10383	Huntleigh Infant Wrap sensor 4 use with Minipulse	\$250.00	EACH
10366	Huntleigh MiniPulse Oximeter with adult probe	\$950.00	EACH
10387	Huntleigh Paeds finger sensor for MiniPulse	\$240.00	EACH
4590	Penlight Disposable	\$45.00	Pkt 6
12230	Pulse Ox- Colour Screen	\$70.00	EACH
9721	Roche Accucheck Performa Meter	\$32.29	EACH
9722	Roche Accucheck Performa Strips	\$44.25	Pkt 50
8977	Roche Accucheck Safe-T-Pro	\$46.46	Box 200
5804	Roche Combur 3 Test E	\$14.50	Pkt 50
5818	Roche Combur 7 Test	\$58.01	Pkt 100
2132	Stethoscope Eartips Soft Black Snap Fit Littman	\$27.00	Set
2067	Stethoscope Littmann Classic II SELong Black Paed	\$196.58	EACH
1054	Stethoscope Nurses Black Dual Head	\$20.00	EACH
40-12-1TP	Welch Allyn Flexiport Cuff-Lrge Adult-1tube-Sphyg	\$41.00	EACH
5041	Welch Allyn Sphygmomanometer D55 - Blue + Cuff	\$170.00	EACH
40509827	Welch Allyn Sphygmomanometer Handheld- Adult Cuff	\$191.00	EACH
DRAPES STERILE			
277547	Foliodrape Sheet 50x50 cm Sterile	\$65.00	Pkt 100
277502	Foliodrape Sheet 75x90 cm Sterile	\$45.00	Pkt 35
2775052	Foliodrape Sheet 95x120 cm Sterile	\$44.00	Pkt 25

CODE	DESCRIPTION	PRICE EXCL GST	UOM
continued. DIAGNOSTICS			
277542	Foliodrape Sheet 90x100 cm Sterile	\$60.00	Pkt 30
EQUIPMENT			
12997	Scale Seca Baby Scales 354	\$240.00	EACH
12993	Seca Carry Bag	\$70.00	EACH
GLOVES			
7035T	Glove Exam Latex Elbow Length Small Sterile	\$9.10	Pair
7036T	Glove Exam Latex Elbow Length Medium Sterile	\$9.10	Pair
7037T	Glove Exam Latex Elbow Length Large Sterile	\$9.10	Pair
1664	Glove Exam USL Latex Powder Free S	\$6.80	Box 100
1665	Glove Exam USL Latex Powder Free M	\$6.80	Box 100
1666	Glove Exam USL Latex Powder Free L	\$6.80	Box 100
1668	Glove Exam USL Latex Powder Free XL	\$6.80	Box 100
1681	Glove Profeel DHD ST P/F 6.5	\$0.95	Pair
1682	Glove Profeel DHD ST P/F 7.0	\$0.95	Pair
1683	Glove Profeel DHD ST P/F 7.5	\$0.95	Pair
1684	Glove Profeel DHD ST P/F 8.0	\$0.95	Pair
1685	Glove Profeel DHD ST P/F 8.5	\$0.95	Pair
19994K	USL Glove Nitrile Polymer Powder Free XS	\$9.90	Box 200
19995K	USL Glove Nitrile Polymer Powder Free S	\$9.90	Box 200
19996K	USL Glove Nitrile Polymer Powder Free M	\$9.90	Box 200
19997K	USL Glove Nitrile Polymer Powder Free L	\$9.90	Box 200
19998K	USL Glove Nitrile Polymer Powder Free XL	\$9.90	Box 200
INCONTINENCE			
16479	MoliCare Premium Lady Pad 1 Drop	\$4.40	Pkt 14
8959	MoliCare Premium Lady Pad 2 Drops	\$6.02	Pkt 14
8960	MoliCare Premium Lady Pad 3 Drops	\$6.96	Pkt 14
1232	MoliCare Premium Lady Pad 4 Drops	\$8.00	Pkt 14
1208	MoliCare Pad 2 Drops	\$7.06	Pkt 28
1209	MoliCare Pad 3 Drops	\$8.42	Pkt 28
1210	MoliCare Pad 4 Drops	\$10.83	Pkt 28
14684	MoliCare Form Normal	\$16.29	Pkt 30
1206	MoliNea Plus 40x 60	\$67.41	CTN 150
1207	MoliNea Plus 60x 90	\$91.08	CTN 100
8204	MoliNea Plus 60x90	\$27.17	Pkt 30
8203	MoliNea Super 40 x 60cm	\$11.60	Pkt 25
14538	MoliCare Premium Fixpants Short Leg Small	\$6.95	Pkt 5
14539	MoliCare Premium Fixpants Short Leg Medium	\$6.95	Pkt 5
14560	MoliCare Premium Fixpants Short Leg Large	\$6.95	Pkt 5
14566	MoliCare Premium Fixpants Short Leg XL	\$6.95	Pkt 5
14536	MoliCare Premium Fixpants Long Leg XXL	\$7.45	Pkt 5
14537	MoliCare Premium Fixpants Long Leg XXXL	\$7.45	Pkt 5
INSTRUMENTS			
600015	Forcep Forester Sponge 24cm Sterile Single Use	\$120.00	Box 20
600016	Needle Holder MayoHegar 14cm Sterile Single Use	\$39.90	Box 20
600018	Scissor Dress Sharp/Blunt 13cm Sterile Single Use	\$22.00	Box 20
529010	Umbilical Scissors 11cm Sterile Single Use	\$35.00	Box 10
100401	Maternity Pack #1 Single Use	\$195.00	Box 10
100402	Maternity Pack #3 Single Use	\$170.00	Box 10
400925	Maternity Scissors 25.5cm Single Use	\$65.00	Box 10
100896	Suture Pack Maternity Single Use	\$175.00	Box 20
WOUNDCARE			
2008	3M Cavilon Durable Barrier Cream 92gm	\$23.54	EACH
8669	3M Cavilon Durable Barrier Cream 28gm	\$8.80	EACH
1286	Dermaplast Sensitive 4cm x 5m	\$19.85	EACH



ORDER BY EMAIL

CUSTOMER SERVICE

customerservices@uslmedical.co.nz

ORDER BY FREEPHONE

AUCKLAND 0800 658 814

CHRISTCHURCH 0800 804 546

USL BIO-MEDICAL SERVICES

AUCKLAND 0800 265 218

EMAIL uslservice@uslmedical.co.nz

AUCKLAND OFFICE

494 Rosebank Road, Avondale



www.uslmedical.co.nz