


# Welch Allyn Green Series™ 900 Procedure Light



6 LEDs

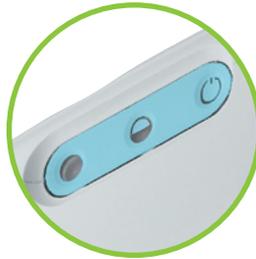

Touchless On/Off  
with Dimming Control

## Whiter, Brighter, and Greener

The GS 900 Procedure Light features:

- Six 5,500° K LEDs that provide a whiter, brighter light than halogen
- Less energy consumption for reduced cost of ownership
- No bulb replacement for cost savings
- Combines outstanding durability, extraordinary maneuverability, unparalleled spot quality and multiple mounting options
- Dimming feature
- Minimal heat dissipation and high color temperature
- 5-year limited warranty

Disposable Sheath  
Available—Part #52630



**WelchAllyn®**  
Advancing Frontline Care™

Green Series Medical Exam Lights



# Welch Allyn Green Series™ 900 Procedure Light

## ORDERING INFORMATION

- 44900 GS 900 Procedure Light with Mobile Stand
- 44900-C GS 900 Procedure Light with Ceiling Mount
- 44900-W GS 900 Procedure Light with Wall Mount
- 52630 Disposable Sheath for GS 600/GS 900 (50/box, 5 boxes/case)

## SPECIFICATIONS

# of LEDs	6 White LEDs
Color Temperature	5,500° K
Light Intensity @ Typical Working Distance	80,000 lux (7,432 footcandles)
LED Life	50,000 hours
Working Distance	24 in (60.9 cm) to 40 in (100 cm)
Spot Size at Working Distance	9 in (23 cm) to 11 in (28 cm)
Electrical	Input: 100-240 V~, 50 Hz-60 Hz, 1.00 A max
Power Cord	#18 AWG, 3-wire grounded cord set, 8' (2.5 m) long
Physical Dimensions	<b>Weight:</b> Total Wall Mount: 17.90 lbs (8.12 kg) Total Mobile Stand: 44.85 lbs (20.34 kg) Total Ceiling Mount: 24.50 lbs (11.11 kg) <b>Luminaire:</b> Height: 13.74 in (34.9 cm) Width: 18.03 in (45.8 cm) Depth: 4.02 in (10.2 cm) <b>Wall to Luminaire Center:</b> 68.9 in (175 cm) <b>Wall Mount Max. Up/Down:</b> 37.4 in (95 cm) <b>Mobile Stand:</b> Height: 68.9 in (175 cm) <b>Mobile Stand Base:</b> Length: 21.97 in (55.8 cm) Width: 22.24 in (56.5 cm) Height: 6.5 in (16.5 cm) <b>Mobile Stand Max. Up/Down:</b> 37.4 in (95 cm) <b>Mobile Stand to Luminaire Center:</b> 37.4 in (95 cm) <b>Ceiling to Luminaire Center:</b> 68.9 in (175 cm) <b>Ceiling Mount Max. Up/Down:</b> 37.4 in (95 cm)
Mounting Options	Wall Mount Mobile Stand Ceiling Mount
Control Box	Touchless On/Off Control with Dimming
Accessory Options	Disposable Sheath #52630
Warranty	5-year limited



## IDEAL FOR USE IN THE FOLLOWING ENVIRONMENTS:

- AMBULATORY CARE CENTERS
- LABOR AND DELIVERY
- HOSPITALS—ER AND ICU
- OUTPATIENT FACILITIES
- OB/GYN
- URGENT CARE CENTERS

## MOUNTING OPTIONS AVAILABLE:



To place an order, contact our Customer Service Department at 1.800.535.6663



**WelchAllyn®**

Advancing Frontline Care™

