



Going further
for health

HARTMANN Folio Surgical Solutions

Providing optimum infection prevention
for your surgery

Foliodrape® Surgical Gowns

FOLIODRESS® GOWN COMFORT

At HARTMANN, we believe that the best surgical gown is one that's actually pleasant to wear, and gives you the confidence of optimum protection against fluid and microbial penetration. That's why we're introducing our range of single-use Foliodress® Surgical Gowns.




Foliodress® Gown Comfort:

- High resistance to fluid penetration
- Comfortable to wear
- Generous cut and length
- Innovative Combitape neck closure for individually adjusted necklines
- Abrasion resistance and low linting
- Sealed seams to reduce the risk of fluid strikethrough
- Anti-static properties
- Strong waist ties and elastic long knit cuffs
- Bi-colour transfer card

Features 4 layers of non-woven polypropylene providing:

- Excellent material stability
- Barrier to bacteria, blood and fluids
- Air permeable / breathable material

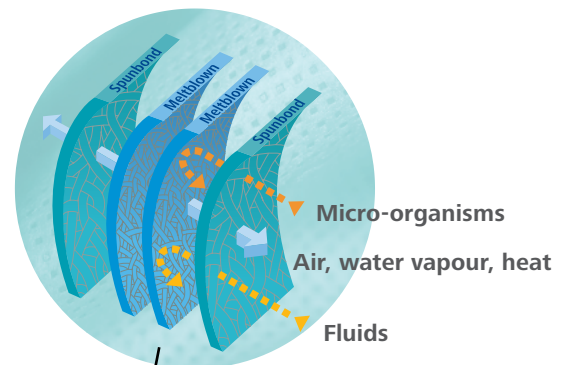


Foliodrape® Surgical Drapes

FOLIODRAPE® DISPOSABLE SURGICAL DRAPES PROVIDE EFFECTIVE DRAPING SOLUTIONS FOR ALL TYPES OF SURGERY

The Foliodrape® Protect and Protect Plus ranges are available in a choice of standard sets, single sterile drapes and supplementary products.

All drapes are low linting, latex free and have anti-static properties. Draping sets are laid out in the order of use for quick and easy handling.

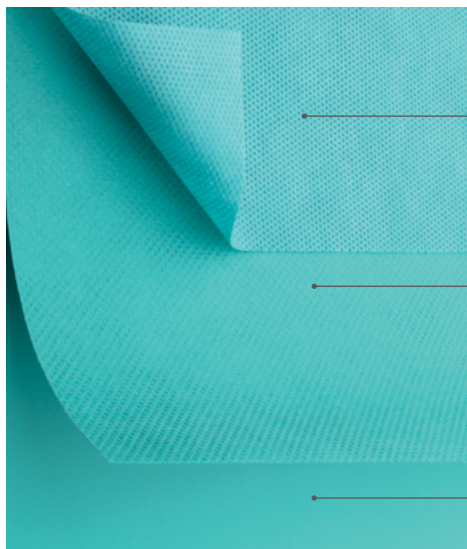

TOP LAYER: HIGHLY ABSORBENT, POLYPROPYLENE REINFORCED ZONE (PROTECT PLUS)

With its optimal size and position on the drape, the reinforced zone ensures reliability and easy handling. It adds material strength to the drape and delivers very high absorption capacity.



MIDDLE LAYER: ABSORBENT POLYPROPYLENE

A special combination of non-woven materials for effective fluid absorption.

BOTTOM LAYER: IMPERMEABLE POLYETHYLENE BARRIER

This layer reduces the risk of post operative infection by preventing the migration of pathogens. Material thickness and weight provides optimal barrier function and tear resistance when wet.

Two layers plus reinforcement zones. Foliodrape® Protect Plus

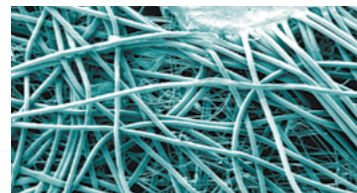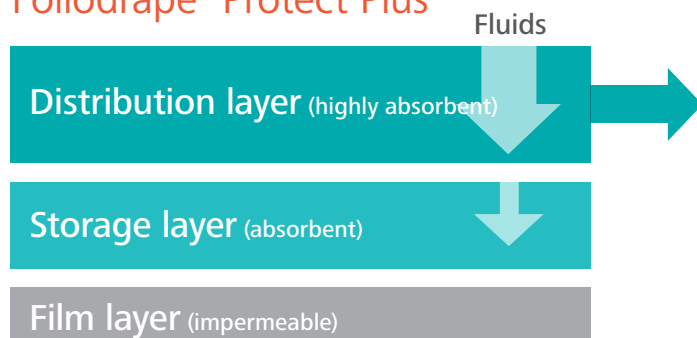


MICROSCOPIC PHOTOGRAPH:

The fluid is quickly absorbed by the material; it is distributed through the fibres in the core (Meltblown) and is safely passed on to the storage layer.

Contents


5	Foliodress® Surgical Gowns	16	Foliodrape® O&G Drapes
5	Foliodress® Comfort Standard Surgical Gown	16	Foliodrape® C-Section Drape
5	Foliodress® Comfort Reinforced Surgical Gown	16	Foliodrape® Lithotomy T Drape
5	Foliodress® Comfort Extra Reinforced Surgical Gown	17	Foliodrape® Under Buttocks Drape
6	Foliodrape® Accessories	17	Foliodrape® Leggings
6	Foliodrape® Absorbent Pad	18	Foliodrape® C-Section Set 3
6	Foliodrape® Tourniquet Cover	18	Foliodrape® Gynae-Abdominal Set
7	Foliodrape® ENT Drapes	19	Foliodrape® Orthopaedic Drapes
7	Foliodrape® ENT Drape	19	Foliodrape® Extremity Drapes
7	Foliodrape® ENT Set	19	Foliodrape® Knee Arthroscopy Drapes
8	Foliodrape® Equipment Covers	20	Foliodrape® Hand Drape
8	Foliodrape® Mayo Stand Cover	20	Foliodrape® Hand Set 1
8	Foliodrape® Protect Back Table Cover	21	Foliodrape® Hand/Foot Set
9	Foliodrape® General Drapes	22	Foliodrape® Hip Drape With Leg Pockets
9	Foliodrape® Surgical Drapes – Non Adhesive	22	Foliodrape® Hip Set 3
9	Foliodrape® Surgical Drapes – Adhesive	23	Foliodrape® Shoulder Arthroscopy Set 1
10	Foliodrape® Surgical Drapes Fenestrated – Non Adhesive	23	Foliodrape® Vertical Isolation Set 1
10	Foliodrape® Surgical Drapes Fenestrated – Adhesive	24	Foliodrape® U Split Drapes
11	Foliodrape® Epidural Drape	24	Stockinette
11	Foliodrape® Protect Plus Drape	25	Foliodrape® Urology
12	Foliodrape® Laparotomy Incise Drape	25	Foliodrape® Urology TUR Drape
12	Foliodrape® Laparoscopy Set 1	25	Foliodrape® Urology Set 1
13	Foliodrape® Universal Set 1	26	Foliodrape® TUR Set 3
13	Foliodrape® Universal Paediatric Set		
14	Foliodrape® Universal Super Set 1		
14	Foliodrape® Universal Set XXL		
15	Foliodrape® Neurology Drapes		
15	Foliodrape® Craniotomy Drape		
15	Foliodrape® Spinal Drape		

Foliodrape® Surgical Gowns

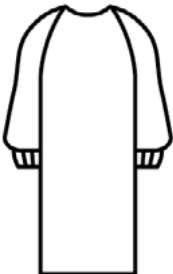	STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	LENGTH	CASE QUANTITY
 Meets EN13795 Standard Performance & AAMI Level 2	927941	81533000	Foliodress® Comfort Standard Surgical Gown	L	130cm	28
	927942	81533010	Foliodress® Comfort Standard Surgical Gown	XL	150cm	28
	927943	81533020	Foliodress® Comfort Standard Surgical Gown	XXL	160cm	28
 Meets EN13795 High Performance & AAMI Level 3	927944	81533030	Foliodress® Comfort Reinforced Surgical Gown	M	120cm	28
	927945	81533040	Foliodress® Comfort Reinforced Surgical Gown	L	130cm	28
	927946	81533050	Foliodress® Comfort Reinforced Surgical Gown	XL	150cm	28
	927947	81533060	Foliodress® Comfort Reinforced Surgical Gown	XXL	160cm	28
 Meets EN13795 High Performance & AAMI Level 3	927948	81533070	Foliodress® Comfort Extra Reinforced Surgical Gown	L	130cm	28
	927949	81533080	Foliodress® Comfort Extra Reinforced Surgical Gown	XL	150cm	28
	927950	81533090	Foliodress® Comfort Extra Reinforced Surgical Gown	XXL	160cm	28


Foliodrape® Accessories

reinforced

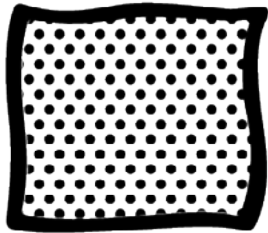

Foliodrape® Absorbent Pad

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
253456	634065	Foliodrape® Absorbent Pad	60 x 75cm	57
258133	636726	Foliodrape® Absorbent Pad	100 x 100cm	45

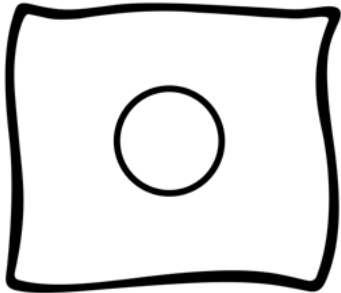


Foliodrape® Tourniquet Cover

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
258505	631311	Foliodrape Tourniquet Cover	45 x 45cm	600

❖ reinforced ❖


Foliodrape® ENT Drapes

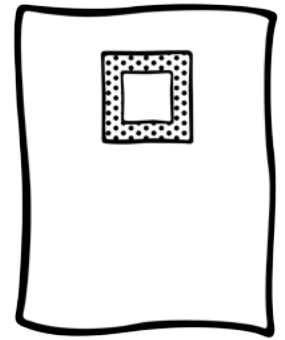

Foliodrape® Protect Plus ENT Drape

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
927901	81532460	Foliodrape Protect Plus ENT Drape	200 x 300cm, Fenestration 20 x 20cm	16

Foliodrape® Protect ENT Set

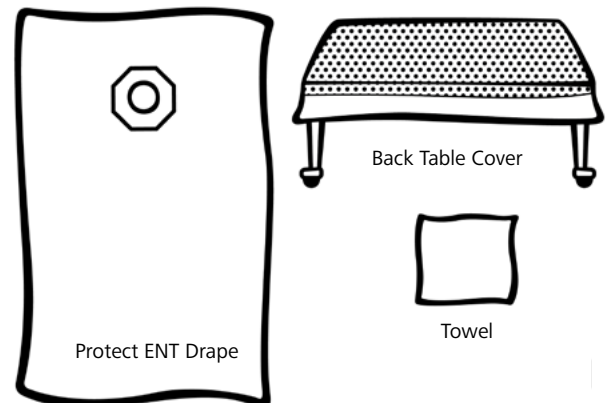

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938912		Foliodrape® Protect ENT Set		14
	–	1 Reinforced Back Table Cover	140 x 160cm	
	–	1 Foliodrape Protect ENT Drape	170 x 250cm Adhesive Fenestration ø 10cm	
	899006	1 Cellulose Towel	30 x 33cm	


Foliodrape® Equipment Covers

reinforced



Foliodrape® Mayo Stand Cover


STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
927914	81532100	Foliodrape Mayo Stand Cover	60 x 145cm	60
938851	635630	Foliodrape Mayo Stand Cover	80 x 145cm	34

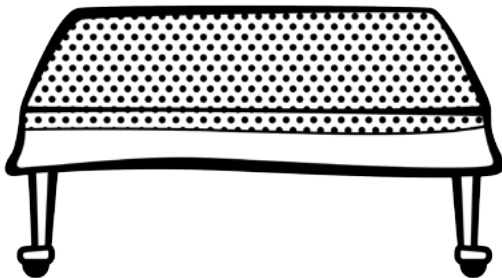

Foliodrape® Back Table Cover

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
927969	81533290	Foliodrape Back Table Cover Reinforced	100 x 150cm	52
250525	631295	Foliodrape Back Table Cover Reinforced	140 x 190cm	24
927905	81532510	Foliodrape Back Table Cover Fully Reinforced	150 x 240cm	18
927906	81532520	Foliodrape Back Table Cover Reinforced	180 x 280cm	24


Foliodrape® General Drapes

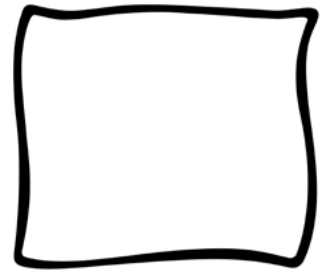

Surgical Drapes - Non Adhesive

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
277500	635865	Foliodrape Protect Drape	45 x 75cm	260
277502	635475	Foliodrape Protect Drape	75 x 90cm	140
277505	637870	Foliodrape Protect Drape	95 x 120cm	100
277506	630200	Foliodrape Protect Drape	90 x 150cm	20
250220	635965	Foliodrape Protect Drape	100 x 150cm	56
927970	81533300	Foliodrape Protect Drape	150 x 100cm	44
277548	635650	Foliodrape Protect Drape	150 x 175cm	36
250224	630480	Foliodrape Protect Drape	150 x 200cm	50
277582	632265	Foliodrape Protect Drape	200 x 170cm	30
277568	638663	Foliodrape Protect Drape	250 x 200cm	17


Surgical Drapes - Adhesive

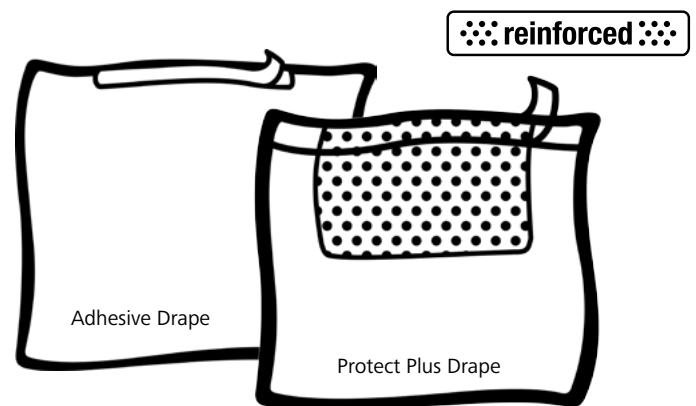

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
277542	632246	Foliodrape Protect Adhesive Drape	90 x 100 cm	120
277543	634761	Foliodrape Protect Adhesive Drape	150 x 175 cm	28
277538	635171	Foliodrape Protect Adhesive Drape	170 x 300 cm	14
938802	635270	Foliodrape Protect Plus Drape	200 x 175 cm	18
938805	637220	Foliodrape Protect Plus Drape	170 x 300 cm	14
927904	81532490	Foliodrape Protect Plus Drape	170 x 300 cm	16


Bold = Adhesive Edge

Foliodrape® General Drapes

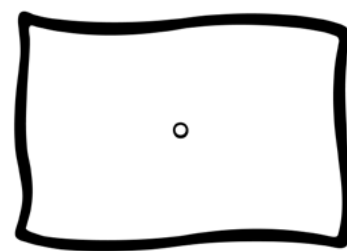

Surgical Drapes - Fenestrated Non Adhesive

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
277503	638380	Foliodrape Protect Fenestrated Drape	45 x 75cm Fenestration ø 7cm	260
277549	632370	Foliodrape Protect Fenestrated Drape	60 x 50cm Fenestration ø 5cm	280
277511	636595	Foliodrape Protect Fenestrated Drape	75 x 90cm Fenestration ø 10cm	160
927968	81533280	Foliodrape Protect Fenestrated Drape	90 x 75cm Fenestration ø 5cm	150

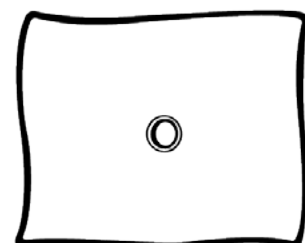


Surgical Drapes - Fenestrated Adhesive

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
277509	637295	Foliodrape Protect Fenestrated Adhesive Drape	45 x 75cm, Fenestration ø variable (2 parts)	260
277510	637290	Foliodrape Protect Fenestrated Adhesive Drape	75 x 90cm, Fenestration ø variable (2 parts)	160
277517	638150	Foliodrape Protect Fenestrated Adhesive Drape	75 x 90cm, Fenestration ø 7cm	160
277555	631960	Foliodrape Protect Fenestrated Adhesive Drape	75 x 90cm, Fenestration ø 10cm	140
277545	634605	Foliodrape Protect Fenestrated Adhesive Drape	120 x 150cm, Fenestration ø 7cm	44
927903	81532480	Foliodrape Protect Plus Fenestrated Adhesive Drape	120 x 150cm, Fenestration □ 15 x 15cm	34

Foliodrape® General Drapes


Foliodrape® Protect Epidural Drape

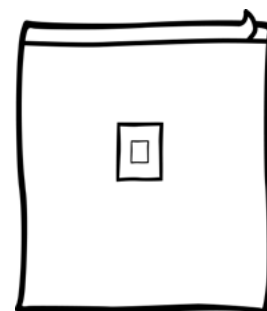

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
927912	81531760	Foliodrape Protect Epidural Drape	125 x 125cm Adhesive Fenestration 10 x 17cm	38
938821	–	Foliodrape Transparent Epidural Drape	90 x 120cm Adhesive Fenestration 10 x 12cm	80

reinforced


Foliodrape® Protect Plus Laparotomy Incise Drape

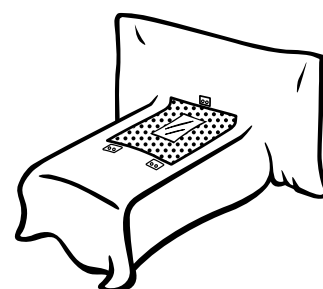

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
927902	81532470	Foliodrape Protect Plus Laparotomy Incise Drape	260 x 200 x 350cm Incise Fenestration 22 x 35cm	8

reinforced


Foliodrape® Protect Plus Laparoscopy Set I

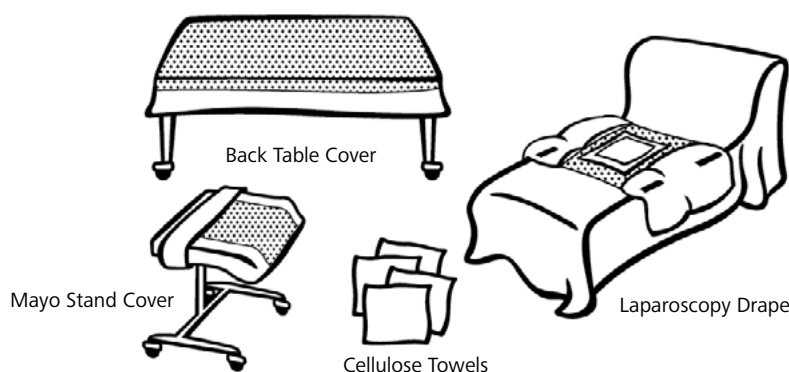

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938709		Foliodrape® Protect Plus Laparoscopy Set I		6
	633965	1 Black Table Cover Reinforced	140 x 190cm	
	630291	1 Mayo Stand Cover	80 x 145cm	
	633457	1 Foliodrape Protect Plus Laparoscopy Drape	250 x 310cm Fenestration 28 x 32cm Bilateral Instrument Pouches	
	899006	4 Cellulose Towels	30 x 33cm	


Foliodrape® General Drapes

Foliodrape® Protect Plus Universal Set I

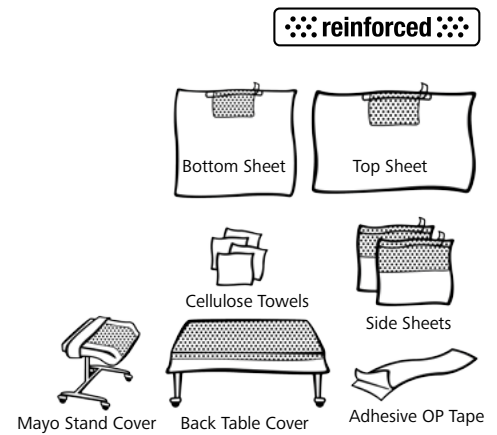

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938712		Foliodrape® Protect Plus Universal Set I		6
	633965	1 Reinforced Back Table Cover	140 x 190cm	
	630291	1 Mayo Stand Cover	80 x 145cm	
	635167	1 Top Sheet	150 x 240 cm	
	635270	1 Bottom Sheet	200 x 175 cm	
	635207	2 Side Sheets	75 x 90 cm	
	997201	1 Adhesive OP Tape	10 x 50cm	
	899006	4 Cellulose Towels	30 x 33cm	

Foliodrape® Protect Plus Universal Paediatric Set

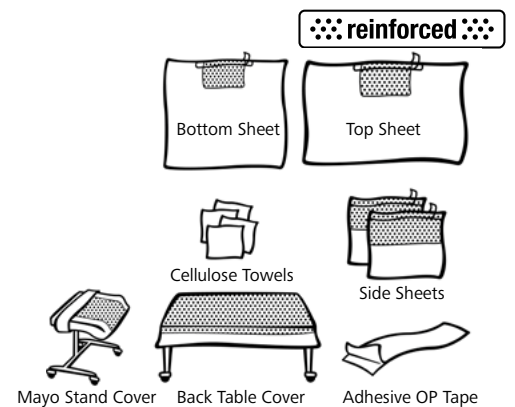


STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938743		Foliodrape® Protect Plus Universal Paediatric Set		6
	633965	1 Reinforced Back Table Cover	140 x 190cm	
	630291	1 Mayo Stand Cover	80 x 145cm	
	435168	1 Top Sheet	150 x 240 cm	
	435271	1 Bottom Sheet	200 x 175 cm	
	435208	2 Side Sheets	75 x 90 cm	
	997201	1 Adhesive OP Tape	10 x 50cm	
	899006	4 Cellulose Towels	30 x 33cm	

Bold = Adhesive Edge


Foliodrape® General Drapes

Foliodrape® Protect Plus Universal Super Set I

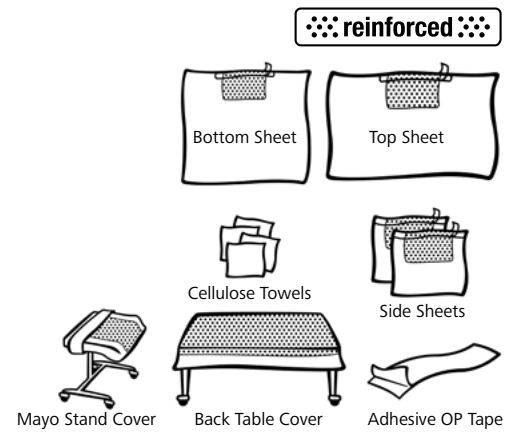

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938707		Foliodrape® Protect Plus Universal Super Set I		5
	633965	1 Reinforced Back Table Cover	140 x 190cm	
	630291	1 Mayo Stand Cover	80 x 145cm	
	637220	1 Top Sheet	170 x 300 cm	
	635270	1 Bottom Sheet	200 x 175 cm	
	635198	2 Side Sheets	90 x 100 cm	
	997201	1 Adhesive OP Tape	10 x 50cm	
	899006	4 Cellulose Towels	30 x 33cm	

Foliodrape® Protect Plus Universal Set XXL

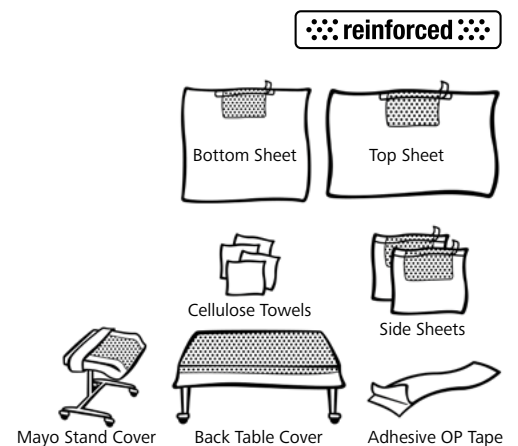


STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938744		Foliodrape® Protect Plus Universal Set XXL		5
	633965	1 Reinforced Back Table Cover	140 x 190cm	
	630291	1 Mayo Stand Cover	80 x 145cm	
	637220	1 Top Sheet	170 x 300 cm	
	634029	1 Bottom Sheet	200 x 200 cm	
	635032	2 Side Sheets	100 x 120 cm	
	997201	1 Adhesive OP Tape	10 x 50cm	
	899006	4 Cellulose Towels	30 x 33cm	

Bold = Adhesive Edge

Foliodrape® Neurology Drapes

reinforced

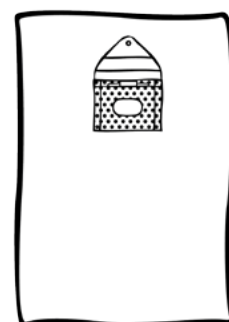

Foliodrape® Protect Plus Craniotomy Drape

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
927966	81533260	Foliodrape Protect Plus Craniotomy Drape	245 x 320cm, Incise Fenestration 40 x 40cm Pouch	8


Foliodrape® Protect Plus Spinal Drape



STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938825	635603	Foliodrape Protect Plus Spinal Drape	320 x 245cm Adhesive Fenestration 9 x 22cm Velcro Strips	8


Foliodrape® O&G Drapes

Foliodrape® Protect C-Section Drape

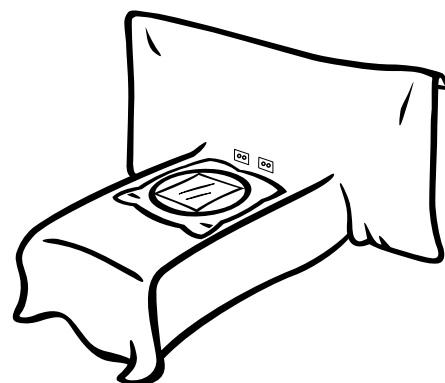

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
927967	81533270	Foliodrape Protect C-Section Drape	260 x 200 x 320cm Full Incise 32 x 32cm Pouch	7

Foliodrape® Protect C-Section Set III

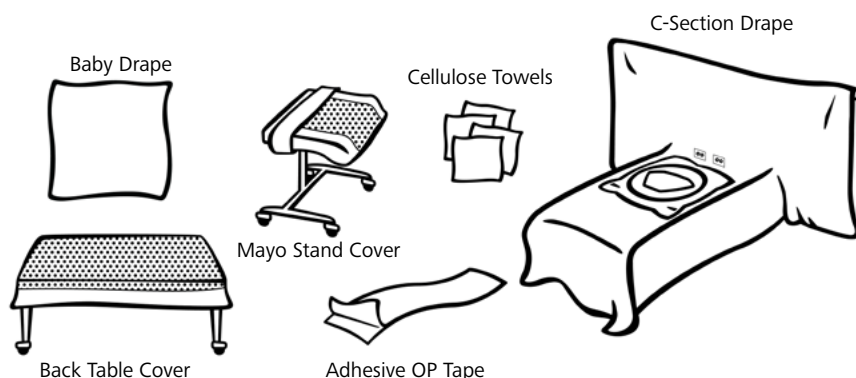


STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938720		Foliodrape® Protect C-Section Set III		5
	633965	1 Reinforced Back Table Cover	140 x 190cm	
	630291	1 Mayo Stand Cover	80 x 145cm	
	333343	1 C-Section Drape	250 x 200 x 315cm Trapezoidal Partial Incise Fenestration 19 x 28cm Pouch	
	—	1 Baby Drape	90 x 100cm	
	997201	1 Adhesive OP Tape	10 x 50cm	
	899006	4 Cellulose Towels	30 x 33cm	

Foliodrape® O&G Drapes

Foliodrape® Protect Plus Lithotomy T Drape




STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
927900	81532440	Foliodrape Protect Plus Lithotomy T Drape	150 x 200cm Fenestration 10 x 15cm	24
927971	81533310	Foliodrape Protect Plus Lithotomy T Drape	200 x 245cm Fenestration 8 x 15cm	18

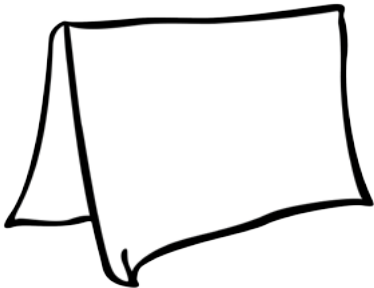

Foliodrape® Protect Leggings

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938842	435084	Foliodrape Protect Legging	75 x 120cm (pack of 1)	20


Foliodrape® Protect Plus Under-Buttocks Drape

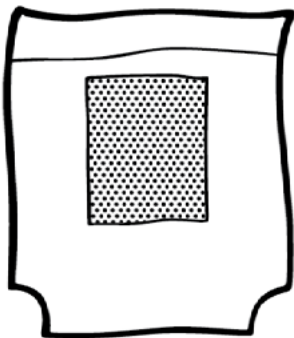

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938829	639936	Foliodrape Protect Plus Under-Buttocks Drape	100 x 120cm	22

Foliodrape® O&G Drapes

Foliodrape® Protect Plus Gynae-Abdominal Set

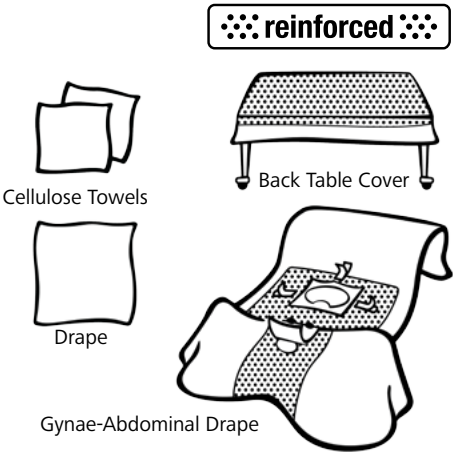


STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938934		Foliodrape® Protect Plus Gynae-Abdominal Set		6
	633965	1Back Table Cover Reinforced	140 x 190cm	
	460039	1 Foliodrape Protect Plus Gynae-Abdominal Drape	230 x 250cm "Kidney Shaped" Fenestration 19 x 29cm Perineal Fenestration 9 x 12cm Integrated Leggings	
	–	1 Drape	90 x 100cm	
	899006	2 Cellulose Towels	30 x 33cm	

Foliodrape® Orthopaedic Drapes


reinforced

Foliodrape® Protect Plus Extremity Drapes

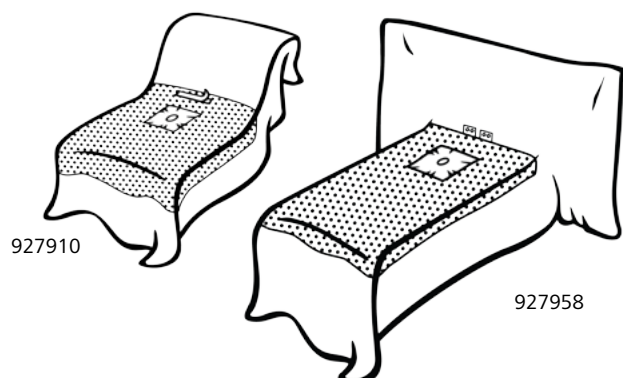

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
927910	81532190	Foliodrape Protect Plus Extremity Drape	218 x 321cm Fenestrated Cuff ø 6cm	7
927958	81533180	Foliodrape Protect Plus Extremity Drape	270 x 225 x 415cm Fenestrated Cuff ø 6cm	7


Foliodrape® Protect Knee Arthroscopy Drapes

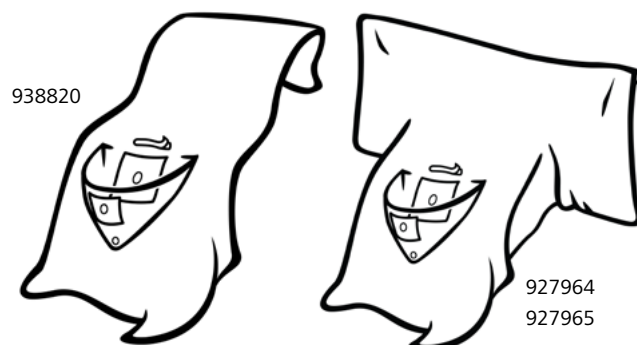

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938820	432803	Foliodrape Protect Arthroscopy Drape	245 x 320cm Fenestrated Cuff ø 6cm Pouch	6
927964	81533240	Foliodrape Protect Arthroscopy Drape	290 x 225 x 335cm Fenestrated Cuff ø 6cm Pouch	6
927965	81533250	Foliodrape Protect Arthroscopy Drape	290 x 225 x 425cm Fenestrated Cuff ø 6cm Pouch	6

Foliodrape® Orthopaedic Drapes


reinforced

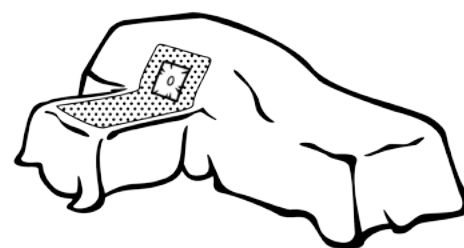

Foliodrape® Protect Plus Hand Drape

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
927962	81533220	Foliodrape Protect Plus Hand Drape	200 x 370cm Fenestrated Cuff ø 3cm	8

reinforced

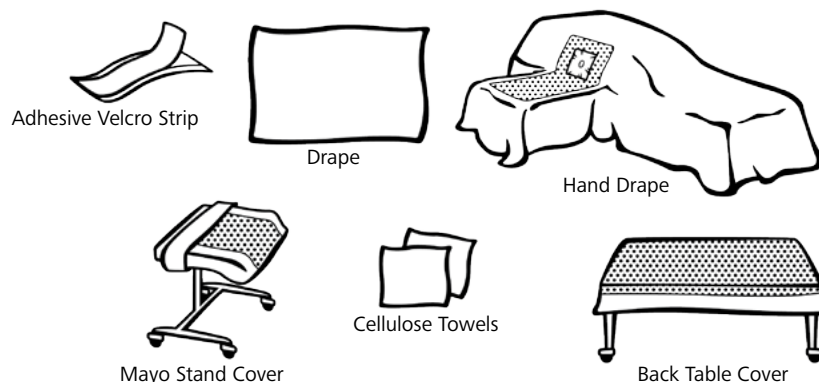


Foliodrape® Protect Plus Hand Set I

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938737		Foliodrape® Protect Plus Hand Set I		5
	633965	1 Back Table Cover Reinforced	140 x 190cm	
	630291	1 Mayo Stand Cover	80 x 145cm	
	—	1 Foliodrape Protect Plus Hand Drape	270 x 370cm Fenestrated Cuff ø 3cm	
	635965	1 Drape	100 x 150cm	
	899006	2 Cellulose Towels	30 x 33cm	
	632383	1 Adhesive Velcro Strip	2 x 33cm	

Foliodrape® Orthopaedic Drapes

reinforced

Foliodrape® Protect Plus Hand/Foot Set

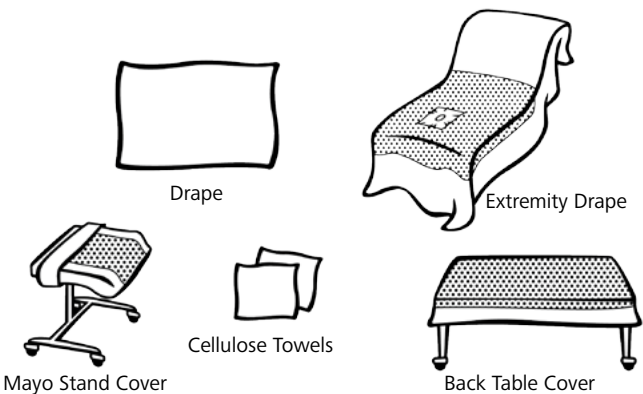


STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938716		Foliodrape® Protect Plus Hand/Foot Set		5
	633965	1 Back Table Cover Reinforced	140 x 190cm	
	630291	1 Mayo Stand Cover	80 x 145cm	
	–	1 Foliodrape Protect Plus Extremity Drape	245 x 320cm Fenestrated Cuff ø 3cm	
	–	1 Drape	100 x 150cm	
	899006	2 Cellulose Towels	30 x 33cm	

reinforced

Foliodrape® Orthopaedic Drapes


Foliodrape® Protect Plus Hip Drape w/leg pockets

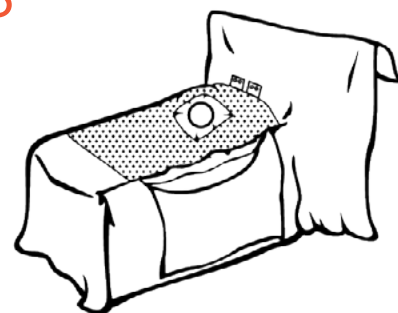

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
927959	81533190	Foliodrape Protect Plus Hip Drape w/Leg Pockets	280 x 225 x 420cm Fenestrated Cuff ø 15cm Bilateral Clear Leg Pockets	6

reinforced

Foliodrape® Protect Plus Hip Set III

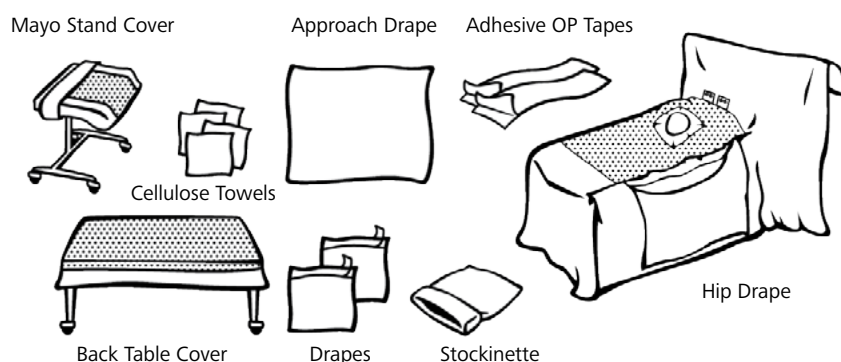


STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938738		Foliodrape® Protect Plus Hip Set III		3
	633965	1 Back Table Cover Reinforced	140 x 190cm	
	630291	1 Mayo Stand Cover	80 x 145cm	
	—	1 Foliodrape Protect Plus Hip Drape	335 x 280cm Fenestrated Cuff ø 6cm Bilateral Clear Leg Pockets	
	637556	2 Drapes	75 x 75cm	
	890744	1 Stockinette	35 x 120cm	
	899006	4 Cellulose Towels	30 x 33cm	
	997201	2 Adhesive OP Tapes	10 x 50cm	
	—	1 Approach Drape	150 x 200cm	

Bold = Adhesive Edge

Foliodrape® Orthopaedic Drapes


reinforced

Foliodrape® Protect Plus Shoulder Arthroscopy Set I

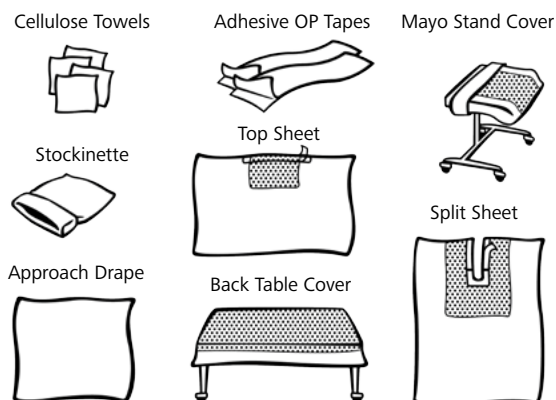

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938702		Foliodrape® Protect Plus Shoulder Arthroscopy Set I		5
	633965	1 Back Table Cover Reinforced	140 x 190cm	
	630291	1 Mayo Stand Cover	80 x 145cm	
	635167	1 Foliodrape Protect Top Sheet	150 x 240 cm	
	333162	1 Foliodrape Protect Split Sheet	200 x 260cm Split 10 x 60cm	
	890745	1 Stockinette	25 x 80cm	
	899006	4 Cellulose Towels	30 x 33cm	
	997201	2 Adhesive OP Tapes	10 x 50cm	
	—	1 Approach Drape	150 x 175cm	


Foliodrape® Vertical Isolation Set I

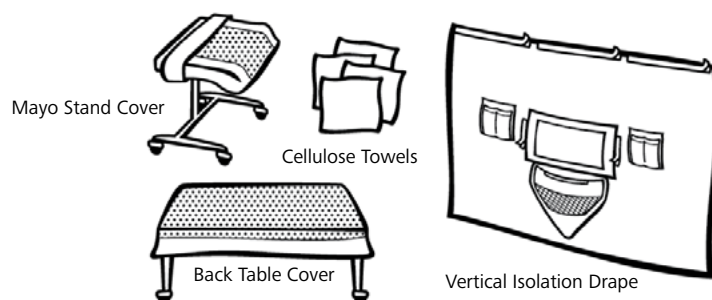

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938705		Foliodrape® Vertical Isolation Set I		6
	633965	1 Reinforced Back Table Cover	140 x 190cm	
	630291	1 Mayo Stand Cover	80 x 145cm	
	333363	1 Foliodrape Vertical Isolation Drape	234 x 320cm Incise Fenestration 30 x 70cm Pouch Instrument Pouches	
	899006	4 Cellulose Towels	30 x 33cm	

Foliodrape® Orthopaedic Drapes


reinforced

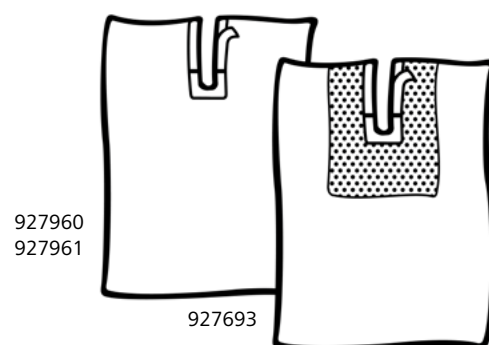

Foliodrape® U and Split Drapes

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
927960	81533200	Foliodrape PE Split Drape	150 x 183cm Split 10 x 52cm	38
927961	81533210	Foliodrape PE Split Drape	150 x 202cm Split 10 x 80cm	35
927963	81533230	Foliodrape Protect Plus U Drape	225 x 300cm Split 15 x 95cm	9

Stockinette

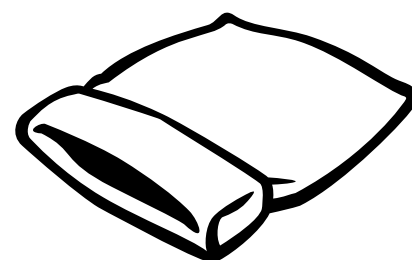

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
927915	81532530	Stockinette Medium	23 x 112cm Knitted Lining Impervious	36
927916	81532540	Stockinette Large	30 x 122cm Knitted Lining Impervious	38


Foliodrape® Urology Drapes

Foliodrape® Protect Urology TUR Drape



STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
927911	81532220	Foliodrape Protect Urology TUR Drape	200 x 260 x 185cm Suprapubic Fenestration 7 x 10cm Genital Fenestration ø7 cm Fingercot Pouch Attached Leggings	10


Foliodrape® Protect Urology Set I

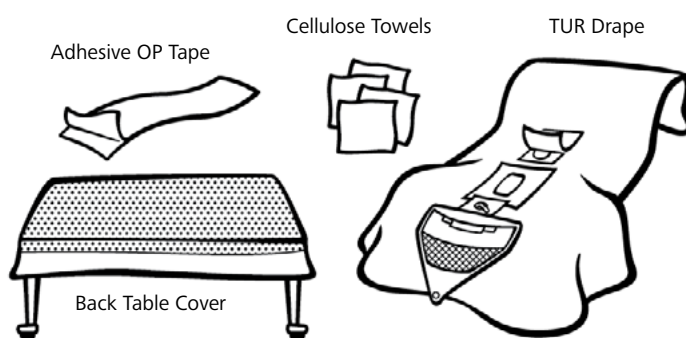

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938700		Foliodrape® Protect Urology Set I		8
	633965	1 Back Table Cover Reinforced	140 x 190cm	
	333352	1 Foliodrape Protect Urology TUR Drape	220 x 230cm Suprapubic Fenestration ø7cm Genital Fenestration 7 x 10cm Fingercot Pouch Attached Leggings	
	997201	1 Adhesive OP Tape	10 x 50cm	
	899006	4 Cellulose Towels	30 x 33cm	

Foliodrape® Urology Drapes

Foliodrape® Protect TUR Set

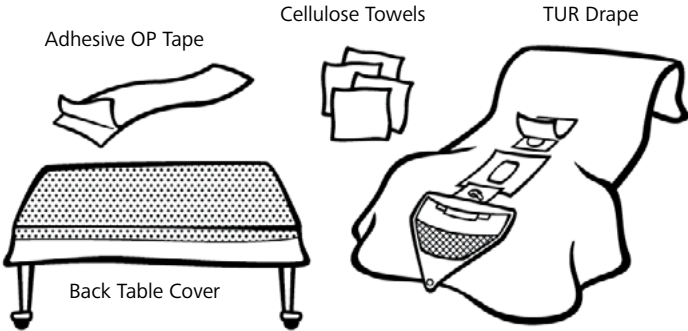

STERILE CODE	NON STERILE CODE	PRODUCT DESCRIPTION	SIZE	CASE QUANTITY
938908		Foliodrape® Protect TUR Set		6
	633965	1 Back Table Cover Reinforced	140 x 190cm	
	–	1 Foliodrape Protect Urology TUR Drape	220 x 230cm Suprapubic Fenestration ø7cm Genital Fenestration 7 x 10cm Fingercot Pouch Attached Leggings	
	997201	1 Adhesive OP Tape	10 x 50cm	
	899006	4 Cellulose Towels	30 x 33cm	

Sterile Code Index

CODE	PAGE
250220	9
250224	9
250525	8
253456	6
258133	6
258505	6
277500	9
277502	9
277503	10
277505	9
277506	9
277509	10
277510	10
277511	10
277517	10
277538	9
277542	9
277543	9
277545	10
277548	9
277549	10
277555	10
277568	9
277582	9
927900	16
927901	7
927902	11
927903	10

CODE	PAGE
927904	9
927905	8
927906	8
927910	18
927911	24
927912	11
927914	8
927915	23
927916	23
927941	5
927942	5
927943	5
927944	5
927945	5
927946	5
927947	5
927948	5
927949	5
927950	5
927958	18
927959	21
927960	23
927961	23
927962	19
927963	23
927964	18
927965	18
927966	14

CODE	PAGE
927967	15
927968	10
927969	8
927970	9
927971	16
938700	24
938702	22
938705	22
938707	13
938709	11
938712	12
938716	20
938720	15
938737	19
938738	21
938743	12
938744	13
938802	9
938805	9
938820	18
938821	11
938825	14
938829	16
938842	16
938851	8
938908	25
938912	7
938934	17



For further information on HARTMANN Surgical products please contact your local Representative or call HARTMANN customer care toll free **1800 805 839**

PAUL HARTMANN Pty Ltd. ABN 35 000 099 589 Level 5, 1 Thomas Holt Drive, Macquarie Park, NSW, 2113.
® Registered trademark PAUL HARTMANN AG.



www.hartmann.com.au