




Rapid  
Simple  
Efficient

 **Haemoband**  
SURGICAL

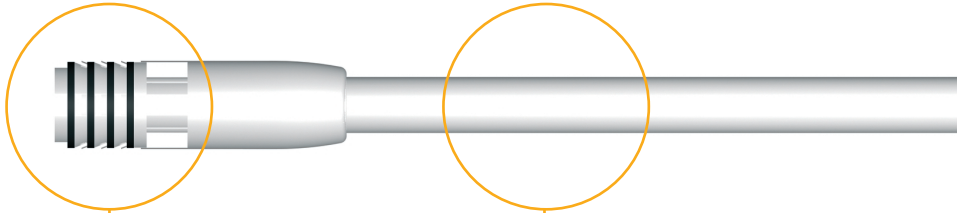

***The multi-ligator tip is pre-loaded with four rubber bands***

Disposable and low cost, the Haemoband offers a very attractive alternative to current devices by providing an easy to use, ergonomically designed product, pre-loaded with four bands.

- Low cost
- Solo operation
- Multi-application
- Rapid
- Accurate


***A slim 120mm insertion probe***

Developed with and by surgeons, the unit has been designed to optimise, simplify and substantially cost-reduce the ligation process.

- Reduced patient discomfort
- Single-use
- Fully disposable unit
- Clinically efficient


ple. Efficient.

The Haemoband Multi-Ligator is a new device that facilitates the easy Rubber Band Ligation (RBL) of haemorrhoids.

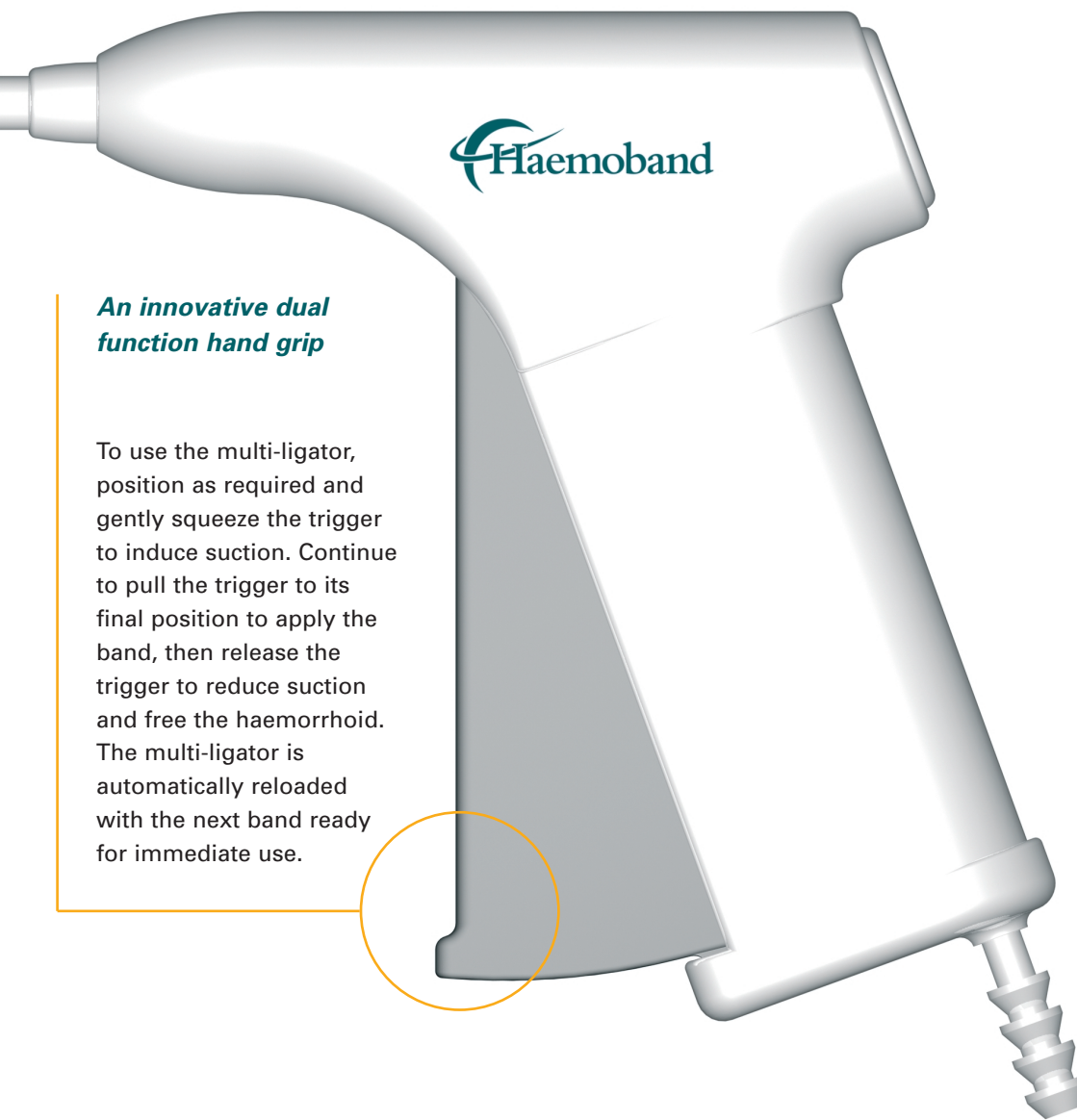

***An innovative dual function hand grip***

To use the multi-ligator, position as required and gently squeeze the trigger to induce suction. Continue to pull the trigger to its final position to apply the band, then release the trigger to reduce suction and free the haemorrhoid. The multi-ligator is automatically reloaded with the next band ready for immediate use.

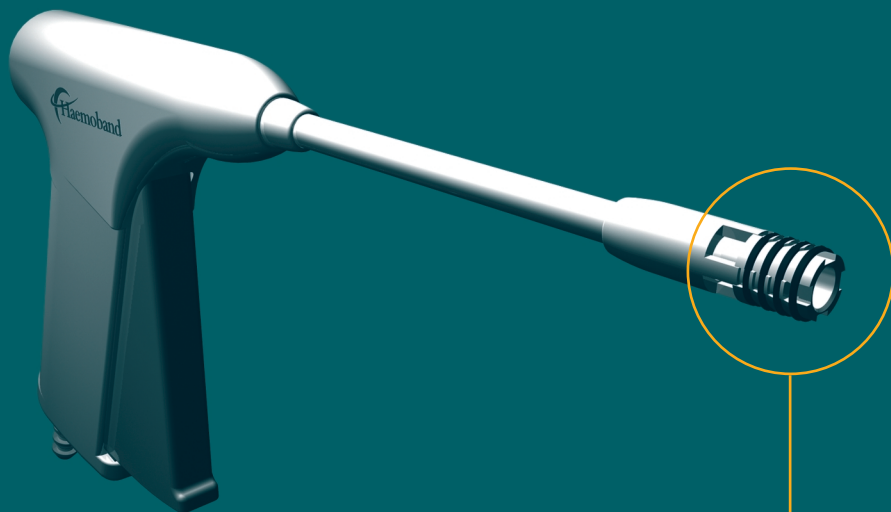

**To be at the forefront in sampling this  
innovative new product for yourself,  
please contact us using the details below:**

Haemoband Surgical Limited  
8 Cromac Quay  
Belfast, BT7 2JD  
Northern Ireland

Tel: +44 (0)28 9033 1200  
Fax: +44 (0)28 9023 9922  
Email: [sales@haemoband.com](mailto:sales@haemoband.com)  
Visit: [www.haemoband.com](http://www.haemoband.com)