



B-lock

F FERNO®

Technical specifications

Weight	0,9 kg	1,9 lbs
Width	300 mm	11,8 inches
Depth	157 mm	6,2 inches
Altezza	240 mm	9,4 inches

Side blocks

Width	85 mm	3,3 inches
Depth	240 mm	9,4 inches
Height	140 mm	5,5 inches

Accessories



Transport and storage bag
Part no. 10-9900-005

 ENTIRELY MADE IN ITALY.



FERNO®

Italy

P +39 051 6860028

F +39 051 6861508

info@ferno.it

www.ferno.it

 800 501711

Schweiz

P +41 41 259 60 00

info@ferno-schweiz.ch

www.ferno-schweiz.ch

FERNO®



B-lock head immobiliser

Lightweight, compact and easy to apply. **Entirely made in Italy**, the B-lock ensures complete and secure immobilisation of the patient's head. Used in combination with the Ferno WizLoc cervical collar, it allows safe and comfortable transportation for the patient.

The B-lock head immobiliser is equipped with two side blocks, a forehead strap and a chin guard, which allow rapid immobilisation thanks to the large Velcro surfaces.

Increased comfort: The material of the B-lock is **soft touch**.



The base of the head immobiliser, in addition to providing stability to the side blocks, has an **integrated QHI** that facilitates head correct positioning.

The B-lock head immobiliser is made of a thermoplastic material with anti-bacterial treatment, which allows easy and optimal cleaning.



The **elliptical ergonomic openings** on the side blocks of the B-lock make it easier to position the head immobiliser.

This makes it possible to maintain direct contact with the patient, who will be more comfortable in the ears' area.




By means of the "H"-shaped Velcro strap, the B-lock can be fixed to the stretcher, thus obtaining **perfect immobilisation** of the patient.



Tested in accordance with the EN 10993 biocompatibility standard.

ERGONOMIC OPENING FOR THE PATIENT'S EARS

INTEGRATED QHI

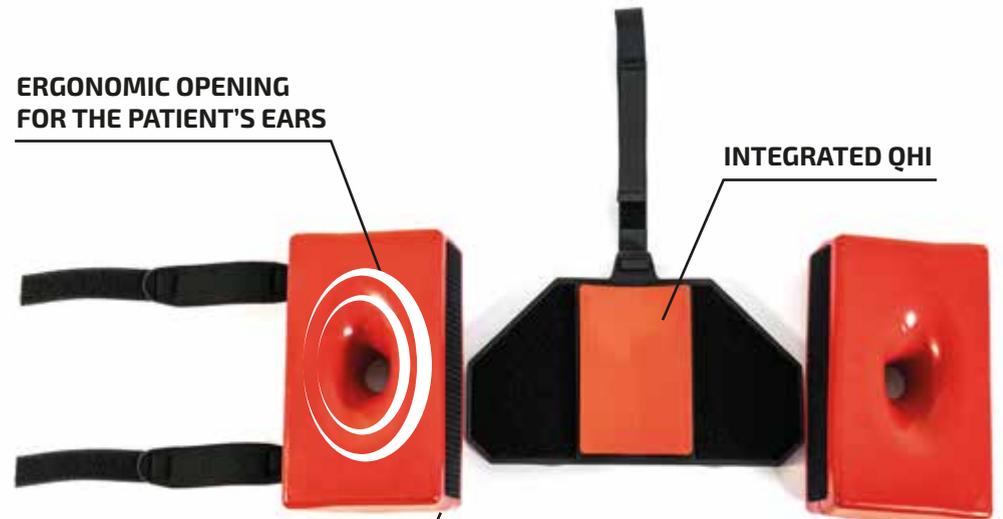

THERMOPLASTIC MATERIAL